TegraDoser-5

Instruction Manual

Manual No.: 15607001

Date of Release: 25.02.2009

TegraDoser-5
Instruction Manual

FCC Notice

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.

Pursuant to Part 15.21 of the FCC Rules, any changes or modifications to this product not expressly approved by Struers A/S could cause harmful radio interference and void the user's authority to operate the equipment.

Always state *Serial No* and *Voltage/frequency* if you have technical questions or when ordering spare parts. You will find the Serial No. and Voltage on the type plate of the machine itself. We may also need the *Date* and *Article No* of the manual. This information is found on the front cover.

The following restrictions should be observed, as violation of the restrictions may cause cancellation of Struers legal obligations:

Instruction Manuals: Struers Instruction Manuals may only be used in connection with Struers equipment covered by the Instruction Manual.

Service Manuals: Struers Service Manuals may only be used by a trained technician authorised by Struers. The Service Manual may only be used in connection with Struers equipment covered by the Service Manual.

Struers assumes no responsibility for errors in the manual text/illustrations. The information in this manual is subject to change without notice. The manual may mention accessories or parts not included in the present version of the equipment.

The contents of this manual are the property of Struers. Reproduction of any part of this manual without the written permission of Struers is not allowed.

All rights reserved. © Struers 2009.

Struers A/S
Pederstrupvej 84
DK-2750 Ballerup
Denmark
Telephone +45 44 600 800
Fax +45 44 600 801

TegraDoser-5

Safety Precaution Sheet

To be read carefully before use

1. The operator(s) should be fully instructed in the use of the machine, any connected equipment and accessories and the applied consumables according to the relevant Instruction Manuals.
2. The machine must be placed on a safe and stable table with an adequate working height. All functions on the machine and any connected equipment must be in working order.
3. Alcohol based consumables: follow the current safety rules for handling, mixing, filling, emptying and disposal of the alcohol-based liquids.
4. If you observe malfunctions or hear unusual noises - stop the machine and call technical service.

The equipment should only be used for its intended purpose and as detailed in the Instruction Manual.

The equipment is designed for use with consumables supplied by Struers. If subjected to misuse, improper installation, alteration, neglect, accident or improper repair, Struers will accept no responsibility for damage(s) to the user or the equipment.

Dismantling of any part of the equipment, during service or repair, should always be performed by a qualified technician (electromechanical, electronic, mechanical, pneumatic, etc.).

Disposal

Equipment marked with a WEEE symbol  contain electrical and electronic components and must not be disposed of as general waste.

Please contact your local authorities for information on the correct method of disposal in accordance with national legislation.

User's Guide

Table of Contents	Page
1. Getting Started	
Checking the Contents of the Packing Box.....	3
Getting Acquainted with TegraDoser-5.....	3
Back of TegraDoser-5.....	4
Unpacking TegraDoser-5.....	5
Placing TegraDoser-5.....	5
Electrical Connections.....	5
Supplying Water (OP versions only).....	5
7-nozzle Dosing Head.....	6
Fitting the Dosing Head to the TegraForce.....	6
Removing the Dosing Head from the TegraForce.....	7
2. Basic Operation of TegraDoser-5	
Front Panel.....	8
Front Panel Controls.....	8
Reading the Display.....	9
Acoustic Signals.....	9
Software Settings.....	10
Changing the Display Contrast.....	10
Editing Numeric Values.....	12
Editing Pre-defined Alphanumeric Values.....	13
Editing New Alphanumeric Values - Using the Text Editor.....	15
Configuring the Suspension and Lubricant Bottles.....	16
Preparation Process Set-up.....	18
Selecting a Preparation Method.....	18
Starting the Process.....	20
Changing parameters while the process is running.....	20
Stopping the Process.....	21
For OP-versions Only.....	21
3. Using Struers Methods	
Introduction.....	22
Viewing Struers Methods.....	22
Viewing a Struers Methods Preparation Method Step.....	23
Copying Struers Methods and Method Steps.....	24

4. Using User Methods	
Introduction	25
Method Groups	25
Creating a New Method Group	26
Renaming Method Groups	27
Deleting Method Groups	28
User Methods	29
Displaying USER METHODS Menu	29
Copying User Methods	30
Inserting/Pasting User Methods	31
Creating a New User Method	32
Renaming User Methods	33
Deleting/Moving User Methods	33
User Method Steps	35
Displaying the User Method Step menu	35
Copying Method Steps	36
Inserting Method Steps	36
Creating a New Method Step	37
Editing Method Steps	38
Defining Surfaces, Lubricants and Suspensions	41
Setting Dosing Levels	43
5. Using Configuration	
Introduction	44
Defining a User Surface Configuration	44
Defining a User Suspension Configuration	46
Defining a User Lubricant Configuration	49
Setting the Operation Mode	51
Changing the Pass Code	54
Resetting the Struers LAN Module	57
6. Trouble-shooting	59
7. Maintenance	
Daily Service	60
Changing Tubes	62
8. Technical Data	65

1. Getting Started

Checking the Contents of the Packing Box

In the packing box you should find the following parts:

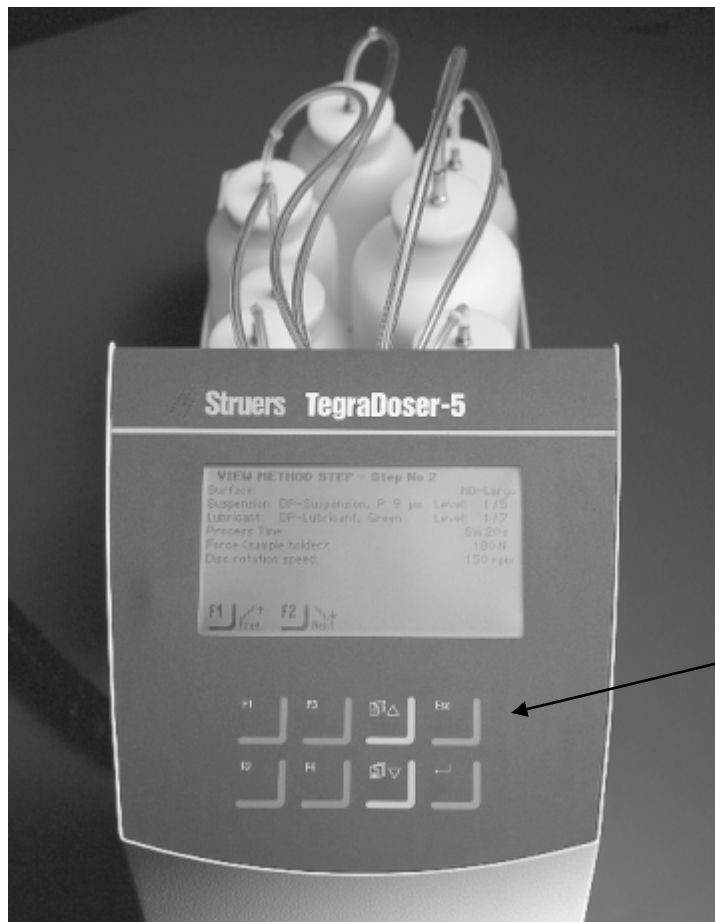
- 1 TegraDoser-5
- 1 Dosing bottle box with bottles and built-in pumps
- 1 Screw
- 2 Silicone tubes for use with alcohol based lubricants
- 1 Set of Instruction Manuals

Note

Check against your order that the machine delivered has the correct number of pumps fitted.
If not, contact your local supplier before using the machine.

Getting Acquainted with TegraDoser-5

Take a moment to familiarise yourself with the location and names of all the TegraPol components:


Front panel control(s)

Back of TegraDoser-5

TegraDoser-5 comes as 2 versions:

- Standard version
- OP-version

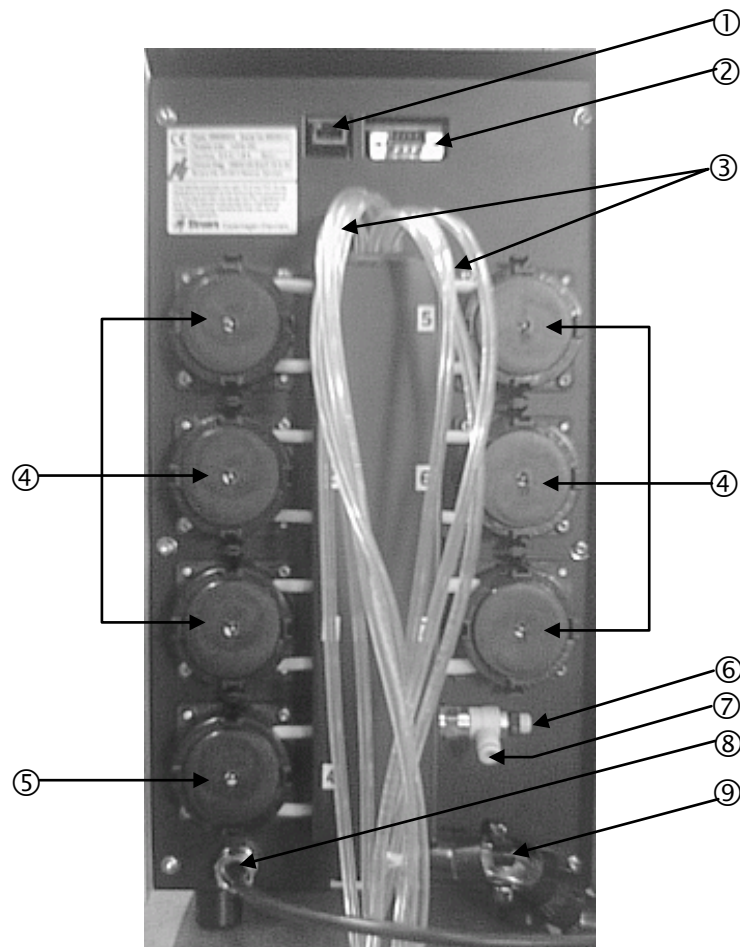Both these versions come as variants, depending on the number of pumps required. Both versions can be fitted with 3 to 6 pumps for diamond suspension or lubricant, only the OP-versions have a pump fitted in position 7 and a connector for the water supply.

- ① Network socket (RJ45 connector)
- ② Service socket (serial port connector)
- ③ Supply tubes from lubricant / suspension bottles
- ④ Pumps (3 to 6 depending on option)
- ⑤ OP-pump (only present on OP-versions)
- ⑥ Regulation valve for water pressure
- ⑦ Water connection from TegraPol for OP-pump (only present on OP-versions)
- ⑧ Power and interface connection from TegraPol
- ⑨ Tubes going to the TegraForce

Unpacking TegraDoser-5

Remove TegraDoser-5 from the foam packing.

Placing TegraDoser-5

Place TegraDoser-5 as close to the TegraPol as possible without it interfering with the operational condition of the TegraPol.

Note

The tubes between TegraDoser-5 and the TegraPol are pre-fitted on delivery. These tubes limit the distance between the two units.

Electrical Connections

There is a single electrical connection between the TegraPol and TegraDoser-5. This provides both the 24V power supply to TegraDoser-5 and a data bus, which enables the two machines to communicate together.

- Switch off the TegraPol.
- Use the single cable attached to the rear of TegraDoser-5 and connect it to the Doser connector on the rear of the TegraPol. To locate this connector, refer to the diagram of the rear of the machine in the TegraPol manual.

**Supplying Water
(OP versions only)**

A water supply is only required when an OP pump is fitted on TegraDoser-5. This comes directly from the TegraPol.

- Push the water tube supplied into the quick fit connector on the rear of the TegraPol.
- Take the other end of the water tube and push it into the quick fit connector on the rear of TegraDoser-5, refer to the diagram of the rear of TegraDoser-5 on page 4.
Do not push the water tube all the way into the connector as that might squeeze the tube and reduce the water flow.
- To adjust the water flow, use the built in regulator at the back of the TegraDoser-5.

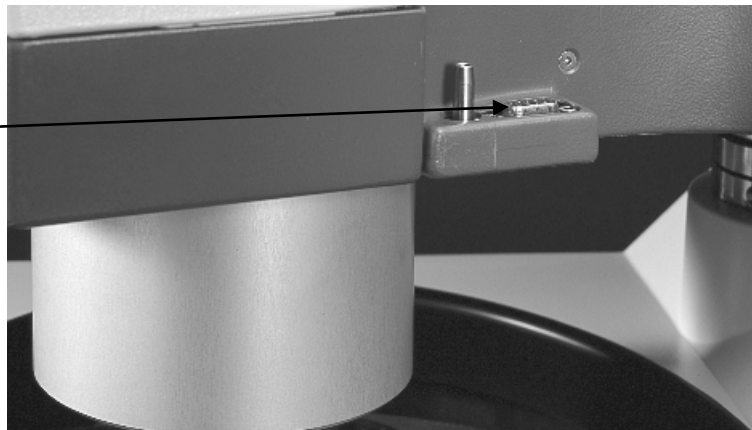
7-nozzle Dosing Head

TegraDoser-5 comes complete with a 7-nozzle dosing head and tubing to allow it to be mounted on the TegraForce.

Note

The length of the tubes between TegraDoser-5 and the dosing head determine the distance between the TegraForce and TegraDoser-5. The factory-mounted tubes are of optimal length and it is not recommended that this length be increased.

D-sub connector
on TegraForce

Fitting the Dosing Head to the TegraForce

- Swing the dosing nozzle arm aside to gain access to the D-sub connector.
- Insert the dosing head into the recess on the side of the TegraForce.
- Press downwards until the dosing head is in position.
- Insert the included screw and tighten it to secure the dosing head.
- Swing the dosing nozzle arm back into position under the body of the dosing head.
- Make sure that the tubes from TegraDoser-5 are not taut so that the TegraForce can be turned freely.

Dosing head
after installation

Dosing nozzle arm

*TegraDoser-5
Instruction Manual*

*Removing the Dosing Head from
the TegraForce*

- Loosen the securing screw located on the top of the dosing head.
- Swing the dosing nozzle arm outwards away from the TegraForce until it is clear of the D-sub connector.
- Grip the top and bottom of the dosing head and pull upwards until it comes free of the TegraForce D-sub connector.
- Ensure that the tubes do not get snagged on the TegraForce as you lift the dosing head clear.




2. Basic Operation of TegraDoser-5

Front Panel

Front Panel Controls

Key	Function
F1	Controls for various purposes. See the bottom line of the individual screens.
F2	Controls for various purposes. See the bottom line of the individual screens.
F3	Controls for various purposes. See the bottom line of the individual screens.
F4	Controls for various purposes. See the bottom line of the individual screens.

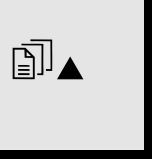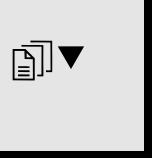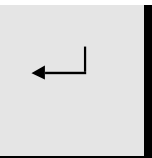Key	Function
	Moves cursor in display upwards or increases the value when editing parameter values.
	Moves cursor in display downwards or decreases the value when editing parameter values.
Esc	Leaves the present menu or aborts functions/changes.
	Enables selected parameter values to be activated for editing. Saves the edited parameter values. Toggles when only 2 parameters available.

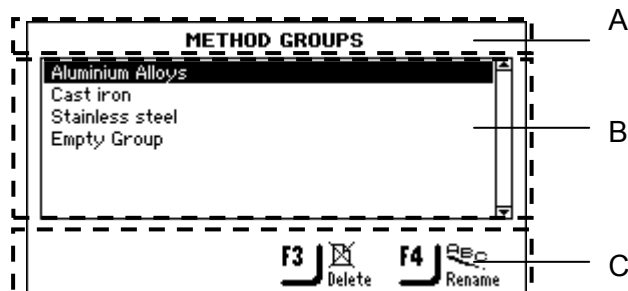
Reading the Display

The display on the front panel provides different levels of status information. For example, when the machine is switched on using the Mains switch located at the rear of the TegraPol, the display informs you which software version is installed:

When operating TegraDoser-5, this display is the user-interface to its software.

The display is primarily divided into 3 areas. The position of these areas and the information they contain are explained in the illustration below, which uses the METHOD GROUPS Menu as an example:

- A** Heading:
this is a navigational aid, telling you where you are in the software's hierarchy.
- B** Information fields:
these will either be numerical values or text fields providing information associated with the process shown in the heading. The inverted text shows the cursor position.
- C** Function key options:
the functions of these keys change depending on the individual situations.

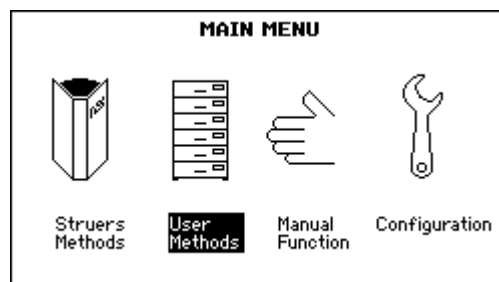
Acoustic Signals

When pressing a key, a short beep indicates that the command has been accepted, whereas a long beep indicates that the key cannot be activated at the moment.

Software Settings

When the machine is started it automatically uses the language selected in the TegraPol configuration menu.

During normal operation, immediately after start up, where the splash screen is displayed, the software goes to the last screen displayed. When starting with the software the MAIN MENU is displayed. This is the highest level in the menu structure. From this menu, you can enter the Struers Methods, User Methods, Manual Function and Configuration menus.

Changing the Display Contrast

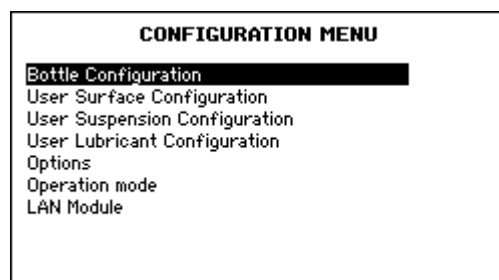
To change the Display contrast setting:

Use the Arrow buttons to select Configuration.


Press ENTER ↵ to activate the CONFIGURATION MENU.


Use the Arrow buttons   to select Options.

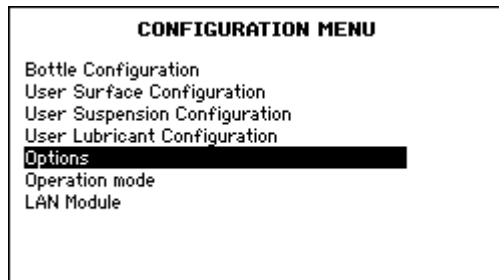TegraDoser-5
Instruction Manual

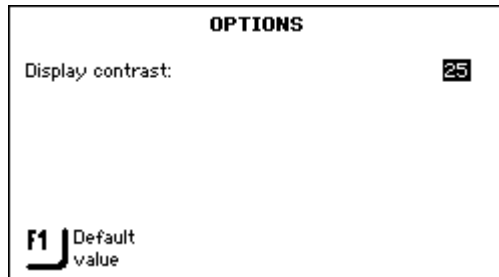Press ENTER ↵ to activate the OPTIONS menu.

Press ENTER ↵ to activate the Display contrast pop-up menu.

Use the Arrow buttons ▲▼ to select the setting you prefer.

Contrast settings are: Default: 25
Contrast range: 0 – 50
Change increment: 1


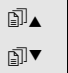
Press ENTER ↵ to accept the setting.

Press Esc 2 times to return to the MAIN MENU.

Editing Numeric Values

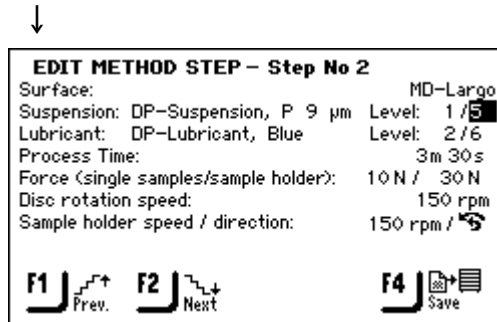 Use the Arrow buttons   to select the value to be changed, e.g. Dosing Level:


↓

 Press ENTER  to edit the value.

↓

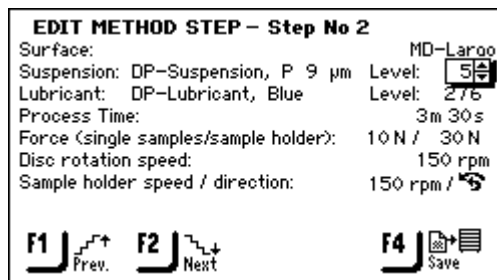A popup box showing the current value and up/down arrows is displayed.



Note:

↓

If there are only two options, the popup box is not displayed. Use ENTER  to toggle between the two options and ignore the following steps.

 Use the Arrow buttons   to increase or decrease the numeric value.

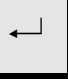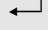
↓

 Press ENTER  to accept the new value.
(Pressing **Esc**, aborts the changes, preserving the original value.)

Editing Pre-defined Alphanumeric Values

If the parameters are already in the system, for example surface types, suspension types and lubricants, then you are presented with a list of possibilities:

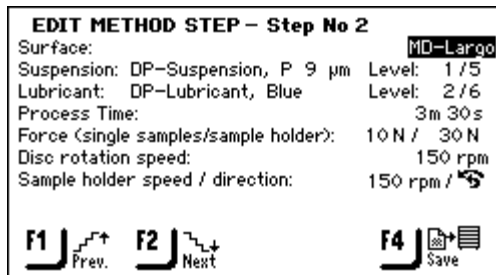 Use the Arrow buttons  ▲ ▼ to select the parameter value to be changed, e.g. Surface type

EDIT METHOD STEP – Step No 2

Surface: MD-Largo

Suspension: DP-Suspension, P 9 µm Level: 1 / 5

Lubricant: DP-Lubricant, Blue Level: 2 / 6

Process Time: 3m 30s

Force (single samples/sample holder): 10 N / 30 N

Disc rotation speed: 150 rpm

Sample holder speed / direction: 150 rpm / ↻

F1 Prev. F2 Next F4 Save

 Press ENTER  to edit the value.

A popup box showing all the options available is displayed:

EDIT METHOD STEP – Step No 2

Surface: MD-Largo

Suspension: Level: 1 / 5

Lubricant: Level: 2 / 6

Process T: 3m 30s

Force (sin: N / 30 N

Disc rotat: 150 rpm

Sample hol: 0 rpm / ↻

SELECT SURFACE TYPE

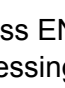
- MD-products
- DP/OP-products
- SP-products
- PG/SiC-paper
- Other products
- User defined

F1 Prev. F2 Next F4 Save

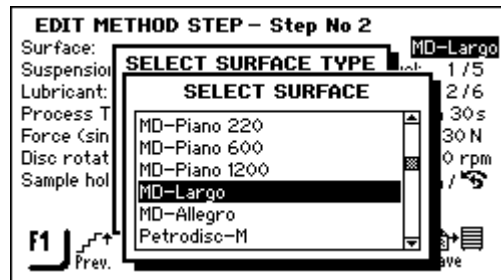
 Use the Arrow buttons  ▲ ▼ to select the correct option.

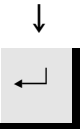
 Press ENTER  to accept the new value.
(Pressing **Esc**, aborts the changes, preserving the original value.)

TegraDoser-5
Instruction Manual

Use the Arrow buttons ▲▼ to select the correct option.

Press ENTER ↵ to accept the new value. (Pressing **Esc**, aborts the changes, preserving the original value.)

Editing New Alphanumeric Values - Using the Text Editor

When new alphanumeric values need to be defined, for example when naming "User defined" consumables (suspensions, lubricants and surfaces) in the CONFIGURATION MENU or renaming Method Groups or User Methods then the TEXT EDITOR is used.

 Use the Arrow buttons  ▲ ▼ to select the value to be changed:

↓

USER SURFACE CONFIGURATION

No.	Surface name	Abr. rule
1	Surf. 1, User type	1
2	Surf. 2, User type	1
3	Surf. 3, User type	1
4	Surf. 4, User type	1
5	Surf. 5, User type	1
6	Surf. 6, User type	1
7	Surf. 7, User type	1

USER METHODS - Aluminium Alloys

AlMg

AlMgSi

AlSi

Empty Method

or

 **F1** Copy

 **F2** Insert

 **F3** Delete

 **F4** Rename

↓

 Press ENTER ↵ to start TEXT EDITOR.

 Press **F4** to start TEXT EDITOR.

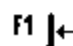
↓


Name field →


Change "Susp 1"
To: Susp 1


TEXT EDITOR

ABCDEFGHIJKLMNOPQRSTUVWXYZ0123456789+*!., :=<>[]'!"?%
 abcdefghijklmnopqrstuvwxyzæøåµ@\
 ÅÅÅÅÇÐÈÉÊËÌÍÎ ÑÒÓÔÕÖÙÚÛÜÝÞß
 ååååçðèéêëìíî ñòóôõöøùúûüýþzef

 **F1** ←

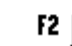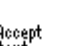 **F2** Accept text

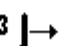 **F3** →

 **F4** Delete

↓

The **F#** key functions are associated with the text displayed in the Name field. Use them to position the cursor in the correct part of the text string, accept and delete the text string displayed.

 Use the Arrow buttons  ▲ ▼ to move around the alphanumeric pad and highlight the correct character.

↓

 Press ENTER ↵ to select the highlighted character.

↓

Repeat the last 2 steps until the desired text string is displayed in the Name field:

 Press **F2** to accept the text string.

Configuring the Suspension and Lubricant Bottles

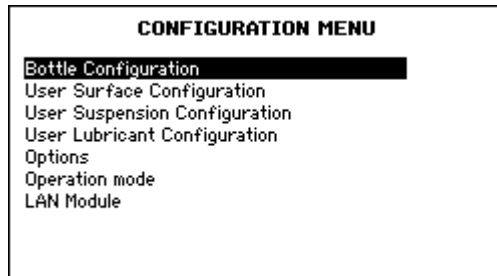
Before TegraDoser-5 can be used in a preparation process, it is necessary to define which suspensions or lubricants are present in the bottles connected to the pumps.

Use the Arrow buttons to select Configuration.

Press ENTER ↵ to activate the CONFIGURATION MENU.

Use the Arrow buttons ⏮ ▲ ▼ to select Bottle Configuration.

Press ENTER ↵ to activate the BOTTLE CONFIGURATION menu.

Bottle	Lub./Susp.	Type
1	Lubricant	???
2	Lubricant	???
3	Lubricant	???
4	Lubricant	???
5	Lubricant	???
6	Lubricant	???
7	Suspension	??? ???

OP-Suspension in bottle 7 only

Use the Arrow buttons ⏮ ▲ ▼ to select the Lub./Susp. column for the individual bottles.


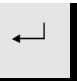
Note!
When present, Bottle 7 can only be configured as OP or AP suspension.

TegraDoser-5
Instruction Manual

 Press ENTER ↵ to toggle between Lubricant and Suspension.




↓
 Use the Arrow buttons ▲▼ to select the Type column.

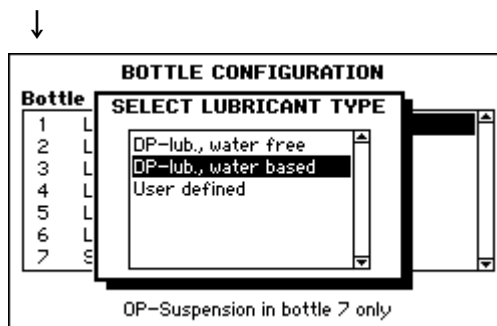↓
 Press ENTER ↵ to display the Lubricant Type or Suspension Type list.

↓
 Use the Arrow buttons ▲▼ to select the correct Lubricant or Suspension Type.

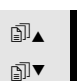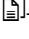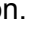↓
 Press ENTER ↵ to display the Lubricants or Suspensions list.

↓
 Use the Arrow buttons ▲▼ to select the correct Lubricant or Suspension.

↓
 Press ENTER ↵ to activate selection.

↓
Repeat the steps above until the contents for all the bottles are defined.

BOTTLE CONFIGURATION			
Bottle	Lib./Susp.	Type	
1	Suspension	DP-P	15 µm
2	Suspension	DP-P	9 µm
3	Suspension	DP-P	6 µm
4	Suspension	DP-P	3 µm
5	Suspension	DP-P	1 µm
6	Lubricant	Green	
7	Suspension	OP-S	

OP-Suspension in bottle 7 only

Press Esc to accept the definitions and return to the CONFIGURATION MENU.

Preparation Process Set-up

Before the TegraPol can be started, a preparation method must be selected. There are two types of methods: Struers Methods and User Methods. Follow the instructions in Chapters 3 and 4, respectively, to select the desired method.

Selecting a Preparation Method

Note!

Depending on the type of TegraForce fitted, the Main Menu on the TegraPol will display different possibilities:
 TegraForce-1: SINGLE SAMPLE PREPARATION is displayed
 TegraForce-3: SAMPLE HOLDER PREPARATION is displayed
 TegraForce-5: SAMPLE HOLDER PREPARATION and SINGLE SAMPLE PREPARATION is displayed

On the TegraPol: if the MAIN MENU is not currently displayed, press **Esc** until you are in the MAIN MENU.

Use the Arrow buttons to select the correct option:

Press ENTER to activate the SINGLE SAMPLE / SAMPLE HOLDER PREPARATION Menu.


SAMPLE HOLDER	
Aluminium Alloys	
Cast iron	
Stainless steel	
Struers Methods	

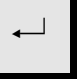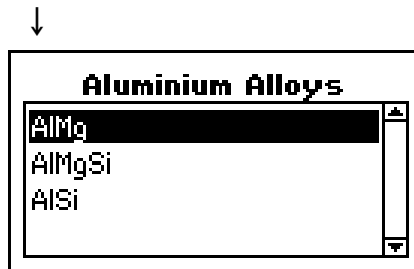
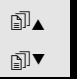
Within this menu, up to 11 groups with preparation methods can be displayed.

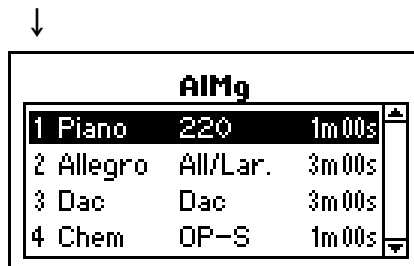
Use the Arrow buttons to select the correct group of methods:

↓
 Press ENTER ↵ to activate the Group menu, which displays a list of up to 20 preparation methods.

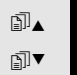
↓
 Use the Arrow buttons   to select the desired preparation method.

↓
 Press ENTER ↵ to activate the Method menu, which can contain up to 20 method steps.

↓ This menu displays the individual steps in the preparation method. The preparation surface, suspension/abrasive and preparation time for each step is displayed.

 Use the Arrow buttons   to select the correct step in the preparation process.

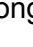
↓
The TegraPol software is now ready to begin the preparation process.

Important!

The TegraForce must be locked in position before the preparation process can be started. An error beep is sounded and an error message is displayed if it is not locked.

Starting the Process

Once the desired preparation method/step has been selected, insert the correct preparation disc on the TegraPol's turntable and press Start  on the TegraPol's control panel to start the turntable rotating.

If water is required to an otherwise dry process, press  to override the process settings and apply water. This occurs as long as  is held.

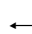
Important!


If you run two successive preparation steps with the same surface, suspension and lubricant, the machine will NOT stop between these two steps, they are considered to be 2 phases of the same step.

However, to stop the machine between two such steps, e.g. to exchange Sic-Paper, select this setting in the TegraPol configuration menu.

Please see the TegraPol Instruction Manual for details.

Changing parameters while the process is running

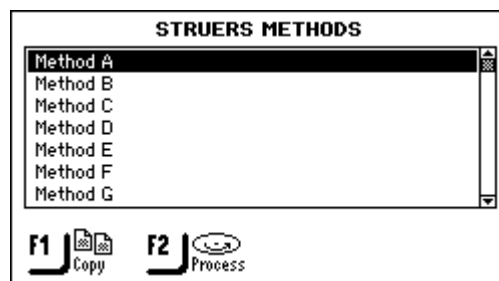While a preparation step is running it might be necessary to change some of the parameters. The only parameters that can be changed in a step that is running are: force and dosing level. To change other parameters the process must be stopped. For fast access to the active method select either Struers Methods or User Methods and press ENTER .

While a method is running, F2 Process provides a shortcut to the step running at the moment. Thus immediate access to the parameters to be changed is possible.

Stopping the Process

To stop the turntable (and water application if using a wet process), press Stop ⏹ on the TegraPol's control panel.

For OP-versions Only

When using an oxide polishing suspension, the last 15 seconds of the preparation time are used to clean the samples and polishing cloth as well as flushing the dosing tubes using water.

If you want to stop this polishing step before it is completed, press Stop ⏹ and the polishing disc and the specimen holder will stop rotating while the water will start/continue to flush for the prescribed 15 seconds. To stop the water flow immediately, press Stop ⏹ again.

Important!

If you stop the OP-version in this manner, ensure that the OP pump is cleaned afterwards to avoid it becoming clogged.

3. Using Struers Methods

Introduction

TegraDoser-5 software is delivered with a set of pre-loaded preparation methods defined by Struers. These methods can be used in their original form or copied and used as a foundation for some user defined methods. The procedures for customising these methods are provided in Chapter 4.

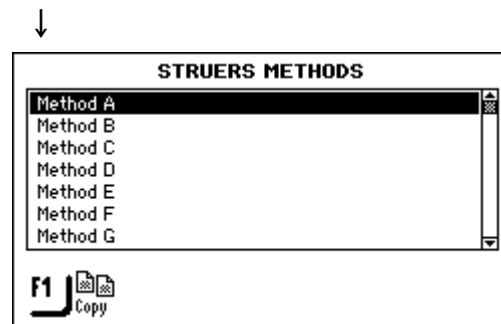
Viewing Struers Methods

Under the Struers Methods menu, the 10 Metalog Guide Methods are listed. To view these methods and the preparation steps they contain:

↓ If the MAIN MENU is not currently displayed, press **Esc** until the MAIN MENU appears.

 Use the Arrow buttons ▲▼ to select STRUERS METHODS:

↓  Press ENTER ↵ to activate the STRUERS METHODS menu.

↓ This menu contains the list of 10 Metalog Guide Methods.

 Use the Arrow buttons ▲▼ to select the correct preparation method:

↓  Press ENTER ↵ to activate the VIEW METHOD - <method name> menu.

↓ This menu displays a list of up to 20 method steps

VIEW METHOD – Method B			
Surface	Susp.	Lub.	Time
1	Primo 220	Water	1 m 00s
2	Largo	DP-P 9 µm Green	5 m 00s
3	Mol	DP-P 3 µm Red	3 m 00s
4	Chem	OP-S	1 m 00s

F1 Copy

Viewing a Struers Methods
Preparation Method Step

↓ With the VIEW METHODS menu for the desired method displayed:

Use the Arrow buttons ▲▼ to select the desired preparation step.

↓

Press ENTER ↵ to activate the VIEW METHOD STEP menu.

↓

VIEW METHOD STEP – Step No 1	
Surface:	MD-Primo 220
Lubricant:	Water
Level:	ON
Process Time:	1 m 30s
Force (single samples/sample holder):	30 N / 180 N
Disc rotation speed:	300 rpm
Sample holder speed / direction:	0 rpm /

F1 Prev. F2 Next

↓

This menu displays the individual preparation parameters for the step: the preparation surface, suspension/lubricant type and settings, downward force applied by TegraForce, turntable rotation speed (and with TegraForce-5 only, direction of samples).

Use the Function buttons **F1** to view the previous step and **F2** to view the next step in the preparation process.

↓

Note

None of the Preparation Methods or Method Steps displayed using the Struers Methods menu option can be edited directly. They must be copied and inserted in the User Methods menu before editing is possible.

Copying Struers Methods and Method Steps

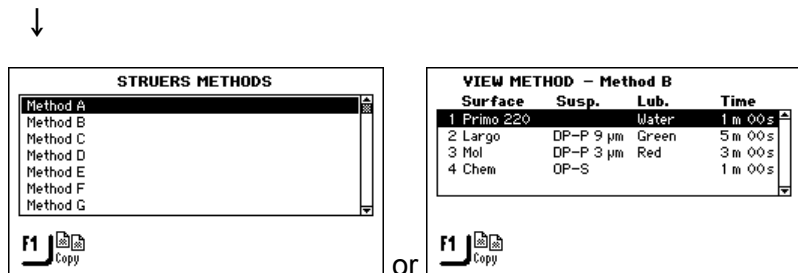
It is possible to copy both the individual methods from the list displayed in the STRUERS METHODS menu and the individual steps displayed in the VIEW METHOD menu.

When a method is copied, all the steps included in this method are transferred across. This provides a complete method that can be used as a foundation and edited to build a new method.

When a step is copied, only the parameter settings for this step are transferred. Steps can only be added to preparation methods that already exist.

To copy Struers Methods and Method Steps:

↓
With the STRUERS METHODS menu or VIEW METHOD menu displayed.
Use the Arrow buttons   to select the desired method or method step.

↓ **Warning!** Using the copy function overwrites whatever is stored on the clipboard.

F1 Use the F1 button to copy the highlighted method or method step.

↓ The copied data is now stored on the clipboard in the random access memory.

Important!

The data on the clipboard will remain there until the copy or delete function is activated again. Using the copy or delete function overwrites whatever exists on the clipboard.

As a method or a method step is stored on the clipboard when deleted this enables you to move a method or a method step from one group to another. Simply delete a method or a single method step and insert it in a different group or preparation method.

4. Using User Methods

Introduction

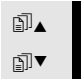
TegraDoser-5 software enables you to build a catalogue of user-defined methods. These can be created either from scratch by defining all the preparation parameters manually, or by copying existing methods which can then be edited via the USER METHODS menu.

To ease the task of managing the preparation methods, they can be grouped together as Method Groups. All the preparation process data available through the User Methods menu can be copied, moved, edited, renamed and deleted.

Method Groups


User Methods that have common attributes, for example they are used on the same class or type of materials, can be grouped together. These are known as Method Groups and are accessed as follows:

↓ If the MAIN MENU is not currently displayed, press **Esc** until the MAIN MENU appears.

↓ Use the Arrow buttons  to select User Methods:

↓ Press ENTER  to activate the METHOD GROUPS menu.

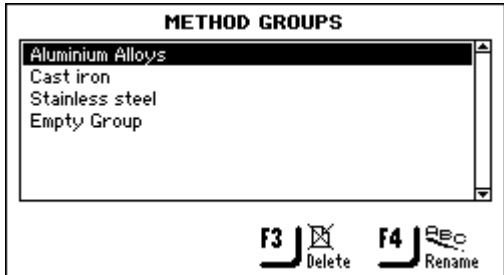↓

The screenshot shows a terminal window titled "METHOD GROUPS". Inside, there is a scrollable list with the following items: "Aluminium Alloys", "Cast iron", "Stainless steel", and "Empty Group". The "Aluminium Alloys" item is highlighted. At the bottom of the window, there are two function key prompts: "F3" with a trash can icon and the word "Delete" below it, and "F4" with a "REC" icon and the word "Rename" below it.

The METHOD GROUPS menu can list up to 10 groups. These are listed alphabetically and from within this menu, you can create new groups, as well as rename and delete existing groups.

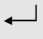
Creating a New Method Group


With the METHOD GROUPS menu displayed.

Use the Arrow buttons ▲▼ to select Empty Group.

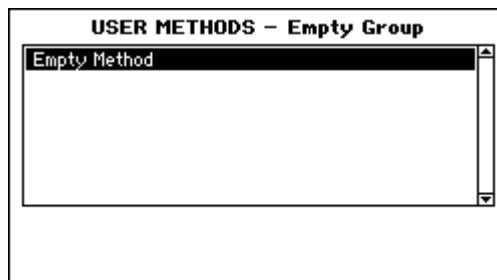Press ENTER  to activate the USER METHODS – Empty Group menu.


In the Empty Group menu, there is only one Empty Method. If you have used the copy function, press F2 to insert the Method currently stored on the clipboard. Otherwise, select the Empty Method for editing (as described in Creating a New User Method on page 32)


As soon as an Empty Method or Group is changed then the name is automatically changed to “Unnamed” Method or Group. The renaming function can then be used to edit the name of the method or group.

Note

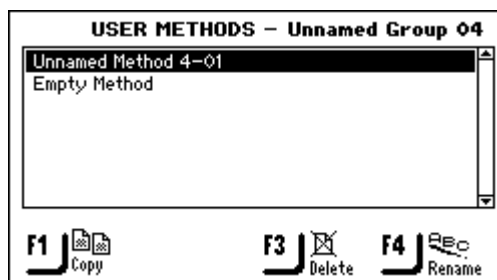At the same time as the Unnamed Method or Group is created, a new Empty Method or Empty Group is created automatically.

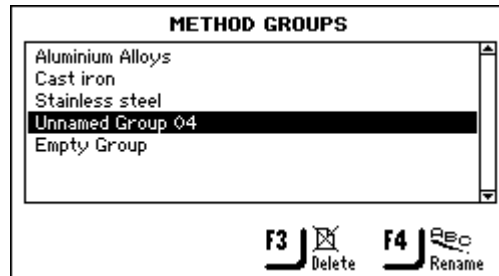
The unnamed group is also included in the METHOD GROUPS menu.

Renaming Method Groups

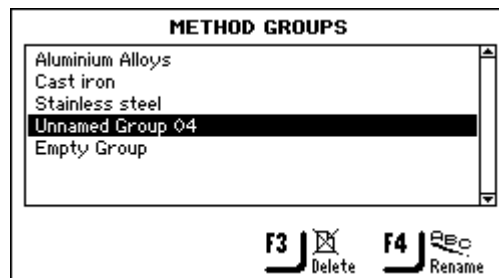With the METHOD GROUPS menu displayed.

New Method Groups have a default new name: e.g. USER GROUPS – Unnamed Group – 04.

Use the Arrow buttons   to select the correct Method Group.

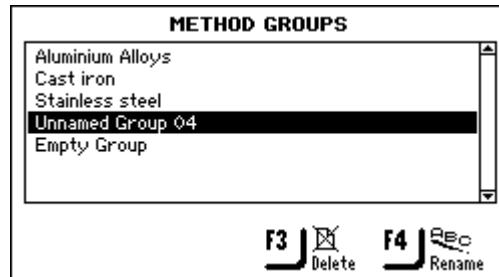
Press **F4** to display the text editor and rename the group.

Deleting Method Groups

Important!
Before a Method Group can be deleted, all the User Methods in this group must be deleted first.

With the METHOD GROUPS menu displayed.

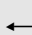
Use the Arrow buttons   to select the correct Method Group.

Press **F3** to delete the selected group.


Before the group can be deleted, confirmation is required
Press ENTER  to confirm the deletion and return to the METHODS GROUP menu.

User Methods

Each Method Group can contain up to 20 User Methods. Unlike the Method Groups, the User Methods are not listed alphabetically. Each User Method provides a series of method steps that are required for the preparation process. Using the USER METHODS menu, you can copy existing User Methods, insert these copied Methods, create new User Methods manually, as well as rename and delete existing User Methods.

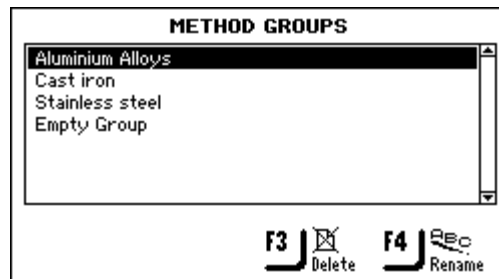
Displaying USER METHODS Menu


With the METHOD GROUPS menu displayed.

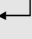Use the Arrow buttons   to select the correct Method Group.

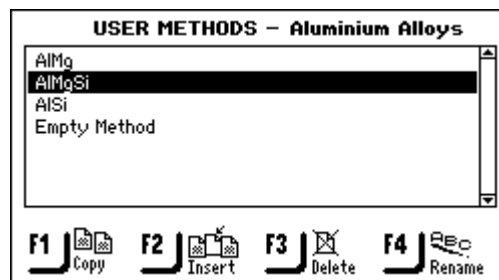Press ENTER  to activate the USER METHODS - <Method group name> menu, e.g. Aluminium Alloys.


Copying User Methods

Like Struers Methods, complete User Methods, including all the process method steps they contain, can be copied and then inserted to create a new User Method.

To copy these User Methods:

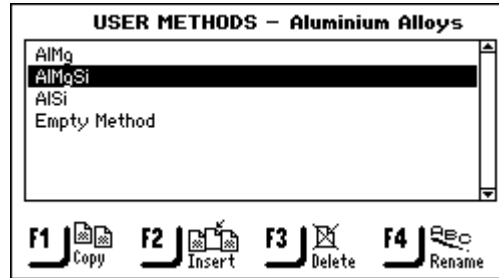With the correct USER METHODS menu displayed.

WARNING!

Using the copy function overwrites whatever is stored on the clipboard.

F1

Use the **F1** button to copy the highlighted method.

The copied data is now stored on the clipboard in the random access memory.

Important!

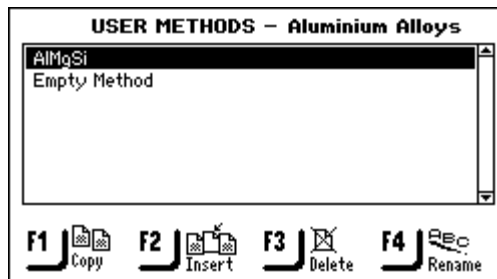
The data on the clipboard will remain there until the copy or delete function is activated again. Using the copy or delete function overwrites whatever exists on the clipboard.

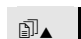
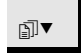
As a method or a method step is stored on the clipboard when deleted this enables you to move a method or a method step from one group to another. Simply delete a method or a single method step and insert it in a different group or preparation method.

Inserting/Pasting User Methods

This feature enables you to insert both Struers Methods and User Methods after they have been copied on to the clipboard. To insert a method:

With the correct USER METHODS menu displayed.

Use the Arrow buttons   to select the position in the list, where the new method is to be inserted.

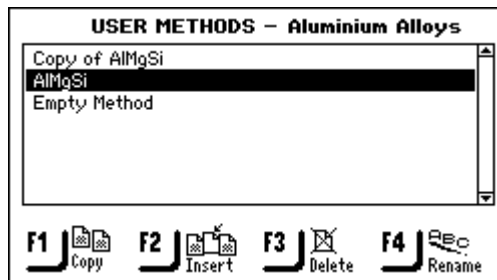
The inserted method will appear above the method highlighted in the list.

Use the **F2** button to insert the contents of the clipboard.

All the method steps associated with this method are also inserted. You are now able to edit the necessary steps – see Editing Method Steps.

Note

It is recommended that you rename the newly inserted method now - use an easy to identify name.

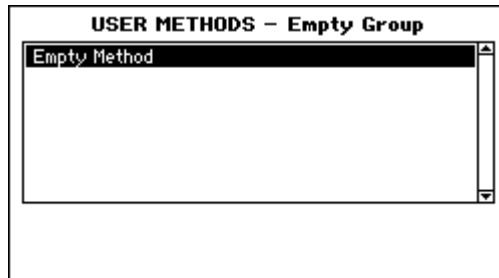
When using the Copy/Insert function the name of the copied method changes automatically to Copy of <method name> when inserted.

If the Delete/Insert function is used the name does not change as the method has been moved, not copied.

Creating a New User Method

To create a new User Method without inserting an existing User Method:

With the correct USER METHODS menu displayed.

↓
When creating a method in a new group, only Empty Method is available.

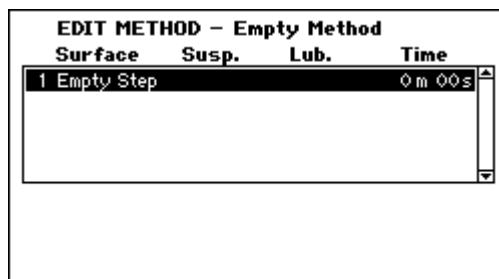
Use the Arrow buttons  to select Empty Method.

↓

Press ENTER  to activate the EDIT METHOD menu.

↓

As soon as this step is changed, the Method name changes to – Unnamed Method - <Group no> 01, and a new Empty Step is created.

If you have used the copy function, press F2 to insert the Method step currently stored on the clipboard.

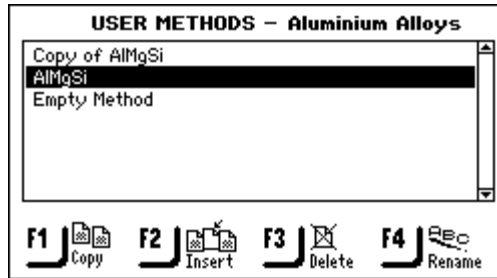
Only a method that contains a step, where the preparation parameters have been set and saved, can be renamed.

Renaming User Methods

Note

A new User Method can only be renamed when it contains at least one method step where the preparation parameters have been set and saved.

With the USER METHODS - < Method group name> menu displayed.

Use New User Methods have a default new name: e.g. Unnamed Method - <Group no> - 01. the Arrow buttons   to select the correct User Method.

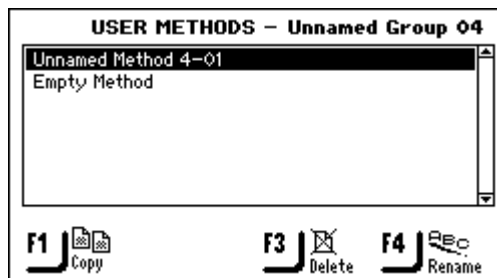
Press **F4** to display the text editor and rename the method.

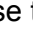
Deleting/Moving User Methods

Important

Once a User Method is deleted it is stored on the clipboard until the next method or method step is deleted or copied. Until then using the Insert function can retrieve it.

With the correct USER METHODS menu displayed.

Use the Arrow buttons   to select the User Method to be deleted.

F3

Use the **F3** button to delete the selected User Method.

All the method steps associated with this method are also deleted.

↓
↵

Press ENTER ↵ to confirm deletion of selected user method.

↓
Esc

Press Esc to go back to the User Method Groups.

↓
▲ ▼

Use the Arrow buttons ▲ ▼ to select the User Method Group where you want to move the method to.

↓
↵

Press ENTER ↵ to activate this User Method Group.

↓
F2

Use the **F2** button to re-insert the deleted User Method.

User Method Steps

Each User Method can contain up to 20 steps. Each Method Step contains a list of process parameters that must be defined and saved before the Method Step can be part of a User Method.

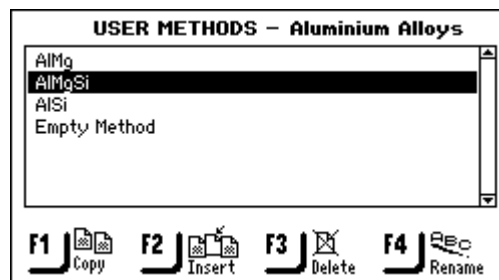
Similar to the User Groups and User Methods, the Copy and Insert functions can be used to allow you to reuse existing Struers and User Method Steps. This means that these are the most fundamental building blocks when creating new User Methods using the insert function.

Important

Two or more method steps can be used as “phases” in a method step. When the same surface, suspension and lubricant are used the other parameters like force and dosing level can be varied to achieve a finer surface finish or to have a more gentle start or finish of the preparation step. These “phases” can have individual preparation times and will be carried out successively without stopping between the “phases”.

Displaying the User Method Step menu


With the correct USER METHODS menu displayed.

Use the Arrow buttons to select the correct User Method.

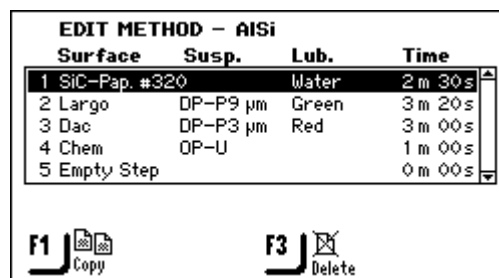Press ENTER to activate the EDIT METHOD - <User Method name> menu.

Copying Method Steps

With the correct EDIT METHOD menu displayed.

EDIT METHOD - AISi			
Surface	Susp.	Lub.	Time
1 SiC-Pap. #320		Water	2 m 30 s
2 Largo	DP-P9 µm	Green	3 m 20 s
3 Dac	DP-P3 µm	Red	3 m 00 s
4 Chem	OP-U		1 m 00 s
5 Empty Step			0 m 00 s

F1 Copy F3 Delete

Use the Arrow buttons ▲ ▼ to select the correct Method Step.

Use the **F1** button to copy the selected User Method Step.

All the parameter settings for this individual step is copied to the clipboard.

Inserting Method Steps

This feature enables you to insert both Struers Method Steps and User Method Steps after they have been copied on to the clipboard. To insert a Method Step:

With the correct EDIT METHOD menu displayed.

EDIT METHOD - AISi			
Surface	Susp.	Lub.	Time
1 SiC-Pap. #320		Water	2 m 30 s
2 Largo	DP-P9 µm	Green	3 m 20 s
3 Dac	DP-P3 µm	Red	3 m 00 s
4 Chem	OP-U		1 m 00 s
5 Empty Step			0 m 00 s

F1 Copy F3 Delete

Use the Arrow buttons ▲ ▼ to select the position in the list, where the new method step is to be inserted.

The new step is inserted above the highlighted step.

Use the **F2** button to insert the contents of the clipboard.

All the parameter settings associated with this method step are also inserted. You are now able to edit the step – see Editing Method Steps.

Creating a New Method Step

With the correct EDIT METHOD menu displayed.

EDIT METHOD – Sintered Carbides			
Surface	Susp.	Lub.	Time
1 Piano 220		Water	1 m 00s
2 Allegro	DiaP. All/Lar.		4 m 00s
3 Dac	DiaP. Dac		4 m 00s
4 Chem	OP-U		1 m 00s
5 Empty Step			0 m 00s

F1 Copy F3 Delete

Use the Arrow buttons ▲▼ to select Empty Step.

Press ENTER ↵ to activate the EDIT METHOD STEP – Step No # (Empty) menu.


EDIT METHOD STEP – Step No 5 (Empty)	
Surface:	MD-Plan
Suspension: DP-Suspension, P 9 µm	Level: 0/0
Lubricant: DP-Lubricant, Blue	Level: 0/0
Process Time:	0m 05s
Force (single samples/sample holder):	10 N / 30 N
Disc rotation speed:	150 rpm
Sample holder speed / direction:	150 rpm /

F1 Prev. F2 Next

To edit the process parameters before pressing F4 and saving the settings, see the next section.

Editing Method Steps

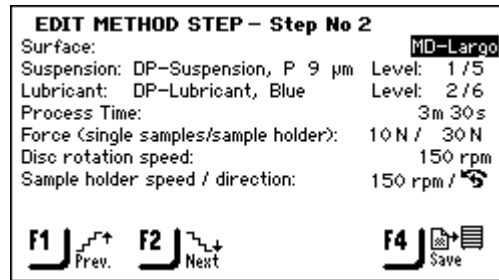Before starting to edit the process parameters for the individual User Method Steps, please note that the parameters displayed vary depending on the type of “Surface:” and “Suspension:” defined. In the screen illustration below, the surface selected requires both an abrasive suspension and a lubricant, but this is not always the situation.


There are 2 types of polishing surfaces: Abrasive Surfaces and Non-abrasive Surfaces. Depending on which one is selected determines what is displayed in the EDIT METHOD menu.

The table below provides a basic overview of the polishing surfaces and the parameters, which are displayed, when they are selected. When “Not required” is indicated in the table, then the line including this parameter, for example “Suspension” is omitted in the screen display.

Type of Surface	Requires	
	Suspension	Lubricant
Abrasive surface	Not required	Water only
Non-abrasive surface	Diamond susp. DiaPro or DiaDuo (All-in-one products)	Not required
	Diamond susp. (water based)	Water-based
	Diamond susp. (water-free)	Water-free
	Oxide susp.	Not required

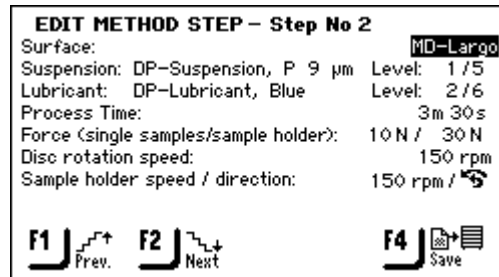With the correct EDIT METHOD STEP menu displayed.

Use the Arrow buttons ▲ ▼ to move between the individual parameters.

Use the **F1** and **F2** buttons to move between the method steps.

When the correct parameter is highlighted.

Press ENTER to edit the parameter.

↓ A popup box showing the current setting and the possible options is displayed. The settings possible for the TegraSystem are illustrated below:

Options	Settings for SINGLE SAMPLES		Change Increment
	TGPol-11, 21 & 31	TGPol-15, 25 & 35	
Speed (rpm)	50 - 300	40 - 600	10 rpm
Process time	60m 00s 05m 00s 01m 00s 00m 05s		30s 10s 5s
Force	10 – 50 for TF-1 only 10 – 100 for TF-5 only		5 N 5 N
Sample holder direction	cw or ccw		
Water	ON or OFF		

Options	Settings for SAMPLE HOLDER		Change Increment
	TegraPol-21 & 31	TegraPol-25 & 35	
Speed (rpm)	50 - 300	40 - 600	10 rpm
Time	60m 00s 05m 00s 01m 00s 00m 05s		30s 10s 5s
Force:	30 - 400		10 N
Dir:	cw or ccw		
Water	ON or OFF		

↓

Use the Arrow buttons  ▲▼ to select the correct value setting or Type of surface, lubricant and suspension.

↓

For Surfaces, Lubricants and Suspension see the following sections for details.

Press the **F4** button to save the edits.

*Defining Surfaces, Lubricants
and Suspensions*

Unlike the other parameter settings, before these parameters can be selected, the type of consumables must be selected first. The same method is used for all three consumable types:

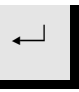
Note


As many different consumables are available, they are sub-divided:
1. Product groups, e.g. MD-Products
2. Product names, e.g. MD-Piano

With the EDIT METHOD STEP menu displayed.


**F1/
F2** Use the **F1** and **F2** buttons to move between the parameters.

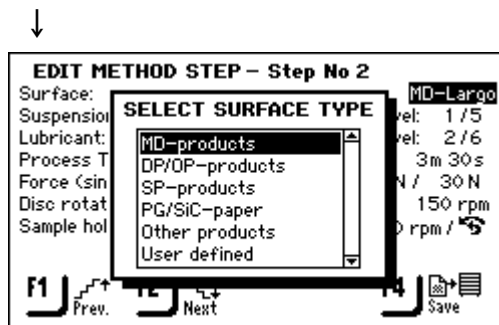↓ When the correct parameter is highlighted.

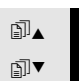 Press ENTER ↵ to activate the SELECT <CONSUMABLE> TYPE menu, e.g. SELECT SURFACE TYPE.


↓

 Use the Arrow buttons   to select the correct <consumable> Type.

↓

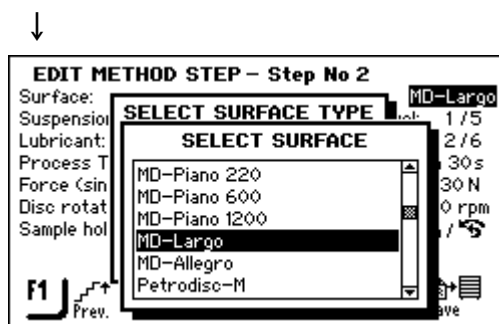
 Press ENTER ↵ to activate the SELECT <consumable> menu, e.g. SELECT SURFACE.

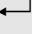Use the Arrow buttons  ▲ ▼ to select the correct consumable.

Press ENTER  to select the highlighted consumable.

Depending on which consumable was selected the EDIT METHOD STEP menu may change its displayed parameters. Use the conventional parameter editing methods to change the remaining parameters.

Setting Dosing Levels

When suspensions and/or lubricants are used in a preparation step, first the type of suspension or lubricant is selected and afterwards the dosing level.

Suspension: DP-Suspension, P 9 µm Level: 1 / 5
Lubricant: DP-Lubricant, Green Level: 2 / 5

Following "Level:", two values can be set, e.g. 1 / 5
 The first value [1] is the pre-dosing level, the amount of suspension or lubricant applied onto the surface before the actual step is started. This is used to provide a lubricated surface to avoid any damage that could occur if the specimens would be running on a dry surface. Depending on the frequency of use and the type of surface different values should be set. For frequently used surfaces a lower value can be used than for surfaces used only once in a while.
 The second value [5] is the dosing level maintained throughout the preparation. This is set according to the type of surface: soft, napped polishing cloths require more lubricant than hard, flat cloths or fine grinding discs. Fine grinding discs require a higher dosing level of abrasive than polishing cloths.

Option	Settings for		Change increment
	Pre-dosing	Dosing	
Dosing Level	0 - 10	0 - 20	1

5. Using Configuration

Introduction

Although the Configuration Menu enables you to change the Display Contrast setting and define the contents of the suspension and lubrication bottles connected to the TegraDoser-5 (both of which are described in Chapter 2), it also allows you to define new user-defined polishing surfaces, lubricants and suspensions.

Defining a User Surface Configuration

Using this option, it is possible to define up to 10 new user-defined surfaces. You are able to define the name, the abrasive rule and lubricant rule for each surface.

To define these parameters:

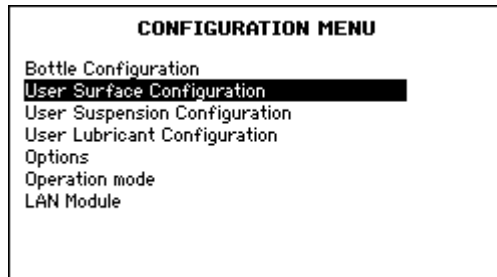
With the MAIN MENU displayed.

Use the Arrow buttons to select Configuration.

Press ENTER ↵ to activate the CONFIGURATION MENU.

Use the Arrow buttons   to select User Surface Configuration.

Press ENTER ↵ to activate the USER SURFACE CONFIGURATION menu.

USER SURFACE CONFIGURATION		
No.	Surface name	Abr. rule
1	Surf. 1, User type	1
2	Surf. 2, User type	1
3	Surf. 3, User type	1
4	Surf. 4, User type	1
5	Surf. 5, User type	1
6	Surf. 6, User type	1
7	Surf. 7, User type	1

Use the Arrow buttons ▲▼ to select the Surface name column opposite the desired No. position.

Press ENTER ↵ to activate the Text editor and define the Surface name.

Use the Arrow buttons ▲▼ to select the Abr. rule column.

Press ENTER ↵ to display the SELECT ABRASIVE RULE menu.

USER SURFACE CONFIGURATION		
No.	Surface name	Abr. rule
1		
2		
3		
4		
5		
6		
7		

SELECT ABRASIVE RULE

1 = Only diamond

2 = Oxide + diamond

3 = Only oxide

4 = Abrasive contained in surface

Use the Arrow buttons ▲▼ to select the correct rule.

Note

It is important to select the correct rules for each surface as the selection will influence the possible choices of suspensions and lubricants available when creating a new preparation method.

If the surface you defined already contains an abrasive, for example SiC paper, select 4 – no abrasives, as this means no abrasives need to be added.

If the surface only is suitable for use with diamond products, select 1.

If both oxide and diamond can be used on the surface, select 2.

If only oxide can be used, select 3.

 Press ENTER ↵ to activate the selected rule and return to the USER SURFACE CONFIGURATION menu.

Defining a User Suspension Configuration

Using this option, it is possible define up to 10 new user-defined suspensions. You are able to define the name, the abrasion type and lubricant rule for each suspension.

To define these parameters:

With the MAIN MENU displayed.

 Use the Arrow buttons to select Configuration.


 Press ENTER ↵ to activate the CONFIGURATION MENU.

CONFIGURATION MENU	
Bottle Configuration	
User Surface Configuration	
User Suspension Configuration	
User Lubricant Configuration	
Options	
Operation mode	
LAN Module	

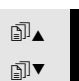 Use the Arrow buttons ▲▼ to select User Suspension Configuration.

 Press ENTER ↵ to activate the USER SUSPENSION CONFIGURATION menu.

USER SUSPENSION CONFIGURATION			
No.	Susp. name	Abr. type	Lub. rule
1	Susp. 1, User type	1	1
2	Susp. 2, User type	1	1
3	Susp. 3, User type	1	1
4	Susp. 4, User type	1	1
5	Susp. 5, User type	1	1
6	Susp. 6, User type	1	1
7	Susp. 7, User type	1	1

Use the Arrow buttons ▲▼ to select the “Susp. name” column opposite the desired No. position.

Press ENTER ↵ to activate the Text editor and define the Suspension name.

Use the Arrow buttons ▲▼ to select the Abr. type column.

Press ENTER ↵ to display the SELECT ABRASIVE TYPE menu.

Note
There are only 2 types of abrasives available.
Be sure to select the correct one.

USER SUSPENSION CONFIGURATION			
No.	Susp. name	Abr. type	Lub. rule
1			
2			
3			
4			
5			
6			
7	Susp. 7, User type	1	1

SELECT ABRASIVE TYPE

1 = Diamond

2 = Oxide

Use the Arrow buttons ▲▼ to select the correct type.

Press ENTER ↵ to activate the selected type and return to the USER SUSPENSION CONFIGURATION menu.

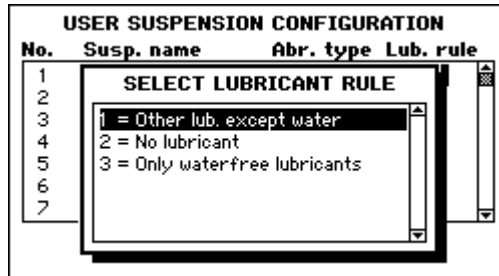
Use the Arrow buttons ▲▼ to select the Lub. rule column.


Press ENTER ↵ to display the SELECT LUBRICANT RULE menu.


Note

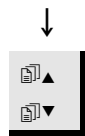If the suspension used requires the use of a lubricant, select 1.
If the suspension doesn't require a lubricant, e.g oxide polishing suspensions or DiaPro, select 2.
If the abrasive is only to be used with water-free lubricants, select 3.

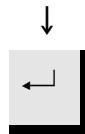Use the Arrow buttons ▲ ▼ to select the correct rule.

Press ENTER to activate the selected rule and return to the USER SUSPENSION CONFIGURATION menu.

Repeat these steps if additional suspensions are to be added. Otherwise follow the next step.

Press Esc to accept the definitions and return to the CONFIGURATION MENU.

Defining a User Lubricant Configuration

Using this option, it is possible to define up to 10 new user-defined lubricants. You are able to define the name and lubricant type for each lubricant.

To define these parameters:

With the MAIN MENU displayed.

Use the Arrow buttons to select Configuration.

Press ENTER ↵ to activate the CONFIGURATION MENU.

CONFIGURATION MENU	
Bottle Configuration	
User Surface Configuration	
User Suspension Configuration	
User Lubricant Configuration	
Options	
Operation mode	
LAN Module	

Use the Arrow buttons ▲▼ to select User Lubricant Configuration.

Press ENTER ↵ to activate the USER LUBRICANT CONFIGURATION menu.


USER LUBRICANT CONFIGURATION		
No.	Lubricant name	Lub. type
1	Lub. 1, User type	2
2	Lub. 2, User type	2
3	Lub. 3, User type	2
4	Lub. 4, User type	2
5	Lub. 5, User type	2
6	Lub. 6, User type	2
7	Lub. 7, User type	2

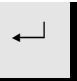
Use the Arrow buttons ▲▼ to select the Lubricant name column opposite the desired No. position.

 Press ENTER ↵ to activate the Text editor and define the Lubricant name.

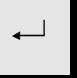↓
 Use the Arrow buttons ▲▼ to select the Lub. type column.

↓
 Press ENTER ↵ to display the SELECT LUBRICANT TYPE menu.

↓

USER LUBRICANT CONFIGURATION		
No.	Lubricant name	Lub. type
1		
2		
3		
4		
5		
6		
7	Lib. 7, User type	2

SELECT LUBRICANT TYPE

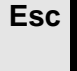
1 = Lubricant containing water
2 = Waterfree lubricant

Note
If the lubricant contains water, select 1.
If the lubricant is water-free and is to be used when preparing water-sensitive materials, select 2.

 Use the Arrow buttons ▲▼ to select the correct type.

↓
 Press ENTER ↵ to activate the selected type and return to the USER LUBRICANT CONFIGURATION menu.

↓
Repeat these steps if additional lubricants are to be added. Otherwise follow the next step.

 Press Esc to accept the definitions and return to the CONFIGURATION MENU.

Setting the Operation Mode

Using this option, it is possible to set various access in the software.
To carry out this operation:

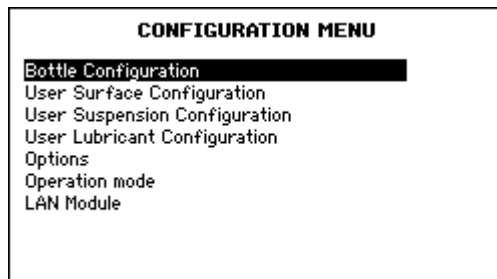
With the MAIN MENU displayed.


Use the Arrow buttons to select Configuration.

Press ENTER ↵ to activate the CONFIGURATION MENU.

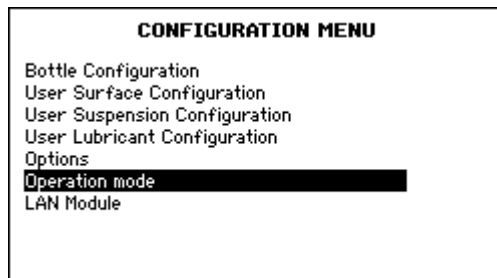Use the Arrow buttons ▲▼ to select Operation mode.

Press ENTER ↵ to activate the Operation mode menu.

TegraDoser-5
Instruction Manual

OPERATION MODE

Current operation mode: Configuration
Pass code: ---

To change the Current operation mode Press ENTER ↵ to enter the Pass code. The default Pass code is 176.

OPERATION MODE

Current operation mode: Configuration
Pass code: 176

Press ENTER ↵ to accept the Pass code.

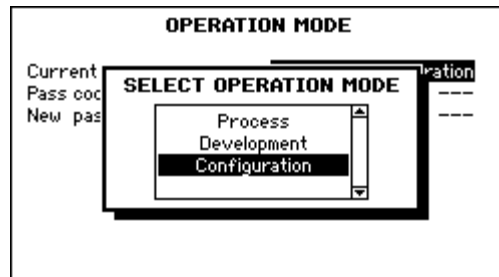
OPERATION MODE

Current operation mode: Configuration
Pass code: ---
New pass code: ---

Press ENTER ↵ to edit the operation mode.

Use the Arrow buttons   to select the Operation Mode.

Operation Modes:

- Process:** Methods can be selected and viewed but no editing is possible.
- Development:** Methods can be selected, viewed and edited
- Configuration:** Methods can be selected, viewed and edited, bottles can be configured and the IP address can be reset.

Changing the Pass Code

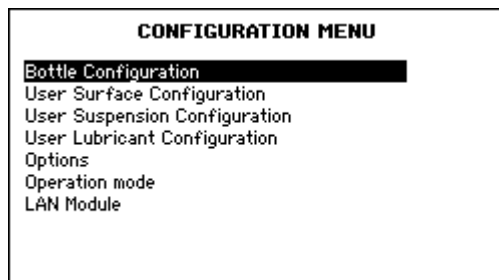
With the MAIN MENU displayed.


Use the Arrow buttons to select Configuration.


Press ENTER ↵ to activate the CONFIGURATION MENU.

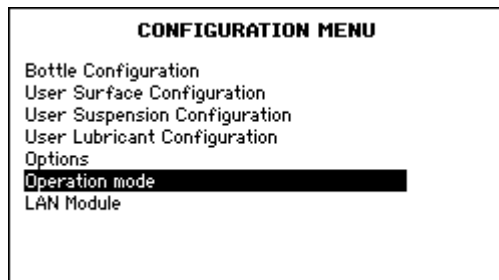Use the Arrow buttons ▲▼ to select Operation mode.

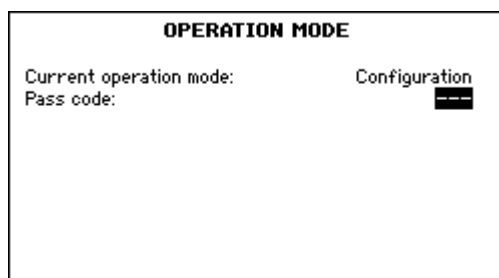Press ENTER ↵ to activate the Operation mode menu.

To change the current pass code Press ENTER ↵ to enter the Pass code. The default Pass code is 176.

OPERATION MODE	
Current operation mode:	Configuration
Pass code:	176

Press ENTER ↵ to accept the Pass code.

Use the Arrow buttons ▲▼ to select: New pass code.

OPERATION MODE	
Current operation mode:	Configuration
Pass code:	---
New pass code:	---

Press ENTER ↵ to access the Pass code.

OPERATION MODE	
Current operation mode:	Configuration
Pass code:	---
New pass code:	0

TegraDoser-5
Instruction Manual

Use the Arrow buttons  ▲▼ to set your personal new Pass code.

Press ENTER  to accept your personal new Pass code.

Resetting the Struers LAN Module

Using this option, it is possible to reset the IP address of the Struers LAN Module to the default address. (Only if the Struers LAN Module is installed).

To carry out this operation:

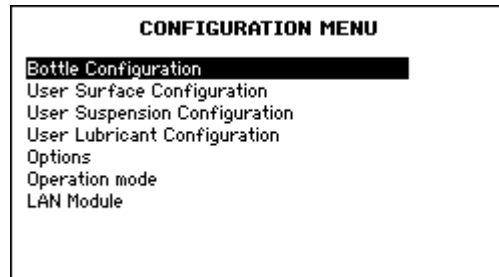
With the MAIN MENU displayed.


Use the Arrow buttons to select Configuration.

Press ENTER ↵ to activate the CONFIGURATION MENU.

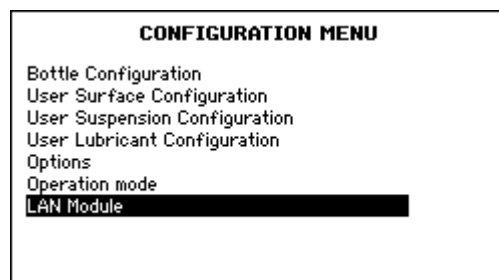
Use the Arrow buttons ↕ to select LAN Module.

Press ENTER ↵ to activate the LAN Module menu.

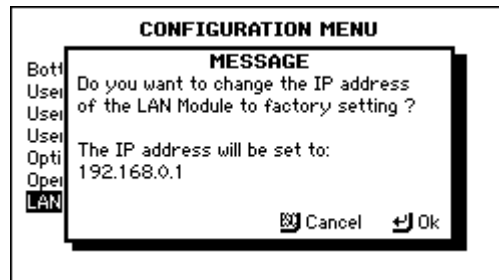*TegraDoser-5
Instruction Manual*

Press ENTER ↵ to reset the IP address to 192.168.0.1

6. Trouble-shooting

Error message	Explanation	Action required
Maximum 20 steps!	Each method can contain a maximum of 20 steps.	
Method database is full!	Database for User Methods is full.	Delete any User Methods that are no longer required.
Method group is full!	Each group can contain a maximum of 20 methods.	Delete an existing method that is no longer required.
Cleaning not allowed during dosing process!	It is not possible to perform cleaning while a process is running .	Clean the tubes after the process is completed.
Suspension and lubricant not compatible!	The selected User Lubricant does not have the correct "Lub. type" or cannot be used together with the selected Suspension. For Example: The selected Suspension cannot be used with water-based Lubricants.	Change the "Lub. type" or choose another Lubricant.
Surface and suspension not compatible!	The selected Suspension does not have the correct "Abr. type" or cannot be used together with the selected Surface. For Example: The selected Surface can only be used with OP- Suspensions and the selected Suspension contains diamonds.	Change the "Abr. type" or choose another Suspension.
Dosing denied, the following suspension is not configured: DP-P 9 µm (Example)	The preparation process cannot start, as the Suspension required is not configured into the TegraDoser-5.	Enter the Bottle Configuration menu and set up the required Suspension. For Example: DP-P 9 µm.
Dosing denied, the following lubricant is not configured: Blue (Example)	The preparation process cannot start, as the Lubricant required is not configured into the TegraDoser-5.	Enter the Bottle Configuration menu and set up the required Lubricant. For Example: Blue.
The method is used for process, certain functions are not allowed!	The attempted function cannot be performed as the method is currently being used for a process For Example: Methods cannot be deleted when they are being used to run a process.	Wait until the process is completed before performing the function.
Process in progress certain functions are not allowed!	The attempted function cannot be performed while a process is running. For Example: Bottle Configuration cannot be changed when a process is running.	Wait until the process is completed before performing the function.

7. Maintenance

Daily Service

Clean all accessible surfaces with a moist cloth.
Check and refill the suspensions and lubricants in the bottles.

Cleaning of Tubes

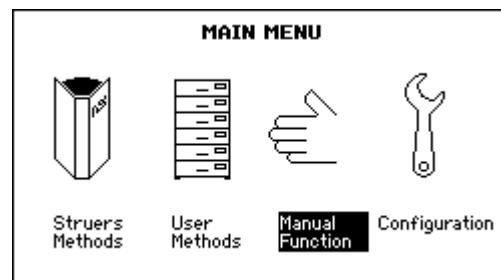
With the MAIN MENU displayed.

Use the Arrow buttons to select Manual Function

Press ENTER ↵ to activate the MANUAL FUNCTION MENU.

CLEANING OF TUBES

No.	Lub./Susp name	Status	Select
1	DiaP. All/Lar.	Used	No
2	DiaP. Largo	Clean	No
3	DiaP. Mol	Clean	No
4	DiaP. Dur	Clean	No
5	DiaP. Dac	Used	No
6	DiaP. Nap-B	Clean	No
7	OP-S	Clean	No

F1 Start cleaning of selected tubes **F2** Select all tubes with Status = Used

The Status column indicates which pumps and tubes have been used

By pressing **F2** the used tubes are selected for cleaning

CLEANING OF TUBES			
No.	Lub./Susp name	Status	Select
1	DiaP. All/Lar.	Used	Yes
2	DiaP. Largo	Clean	No
3	DiaP. Mol	Clean	No
4	DiaP. Dur	Clean	No
5	DiaP. Dac	Used	Yes
6	DiaP. Nap-B	Clean	No
7	OP-S	Clean	No

F1 Start cleaning of selected tubes **F2** Select all tubes with Status = Used

↓ The Select column now shows which pumps are selected for cleaning. Pressing Enter on a highlighted parameter will toggle the selection between Yes and No. Thus you can select/deselect pumps from being cleaned.

F1

Pressing **F1** starts the cleaning sequence

Then simply follow the messages on the screen to clean the pumps.

When the cleaning is completed the selected tubes have their status changed back to Clean.

Now you can also change to a different type of suspension if necessary.


CLEANING OF TUBES			
No.	Lub./Susp name	Status	Select
1	DiaP. All/Lar.	Clean	No
2	DiaP. Largo	Clean	No
3	DiaP. Mol	Clean	No
4	DiaP. Dur	Clean	No
5	DiaP. Dac	Clean	No
6	DiaP. Nap-B	Clean	No
7	OP-S	Clean	No

F1 Start cleaning of selected tubes **F2** Select all tubes with Status = Used

Changing Tubes


When working with alcohol based lubricants, the Novoprene tubes mounted in the pumps supplied with TegraDoser-5 will harden over time. Therefore a piece of silicone tubing is supplied with the TegraDoser-5 as Silicone has a better resistance against alcohol.

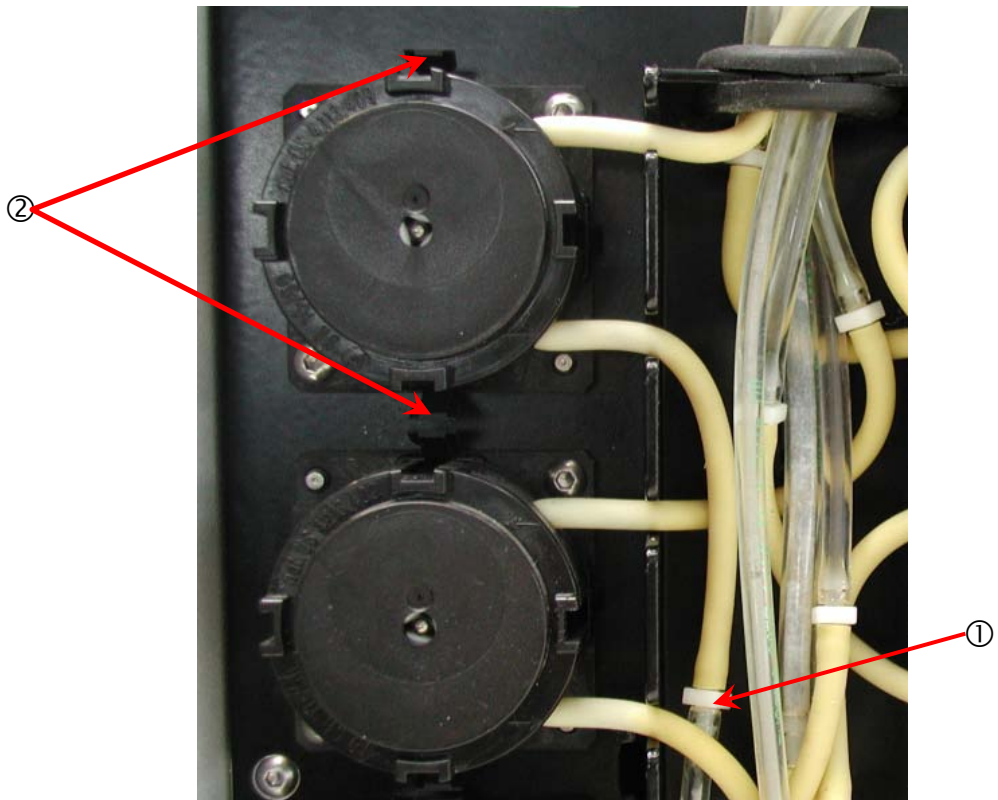To exchange the Novoprene tube with a Silicone tube:

- Separate the doser tubes at the white couplings ① (the couplings should stay on the tubes connected to TegraDoser-5)
- Press the two tabs ② and remove the pump from the axle.

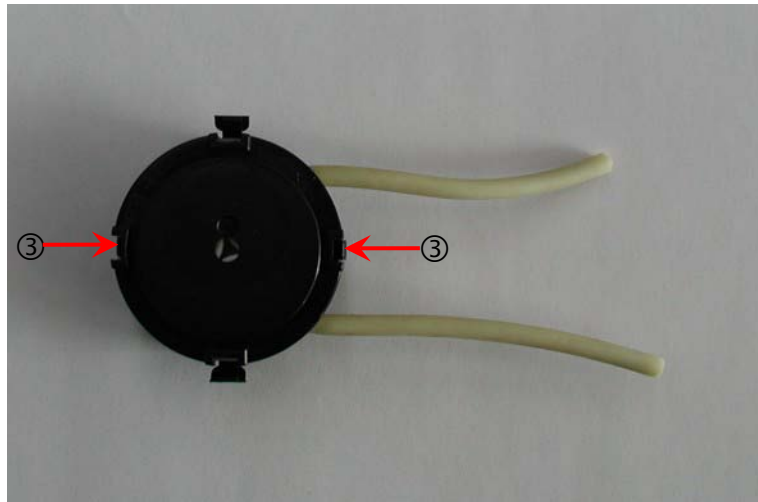- Press the two tabs on the pump ③ and remove the bottom cover.

- Remove the 3 loose rollers and replace the Novoprene tube with the silicon tube.

- Replace the 3 rollers in the pump housing.

- Re-mount the bottom cover.
- Re-connect the tubes to the tubes on the TegraDoser-5 and press the pump back onto the axle.
- Check that the tubes are connected correctly so that liquid is pumped to the head.

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Date of Release 01.06.2009
Update 11.08.2010

TegraDoser-5
Spare Parts and Diagrams

**Always state *Serial No* and *Voltage/frequency*
if you have technical questions or when ordering spare parts.**

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


Some of the drawings may contain position numbers not used in connection with this manual.

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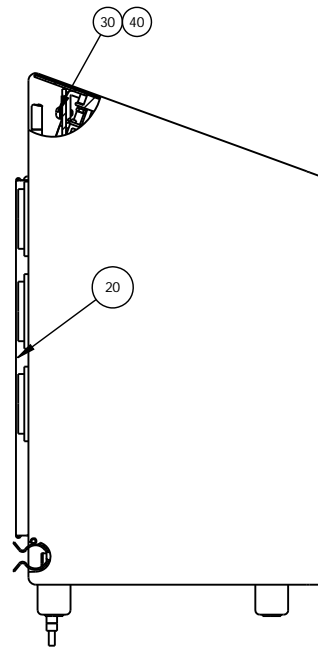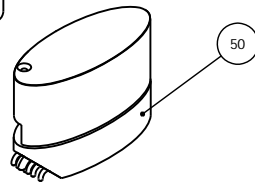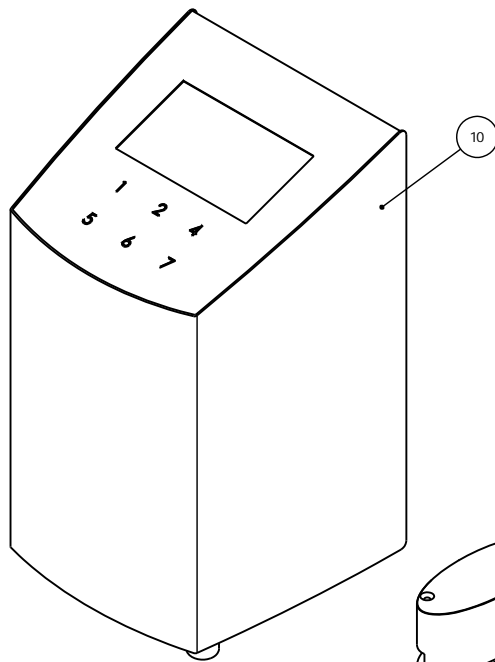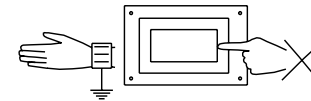
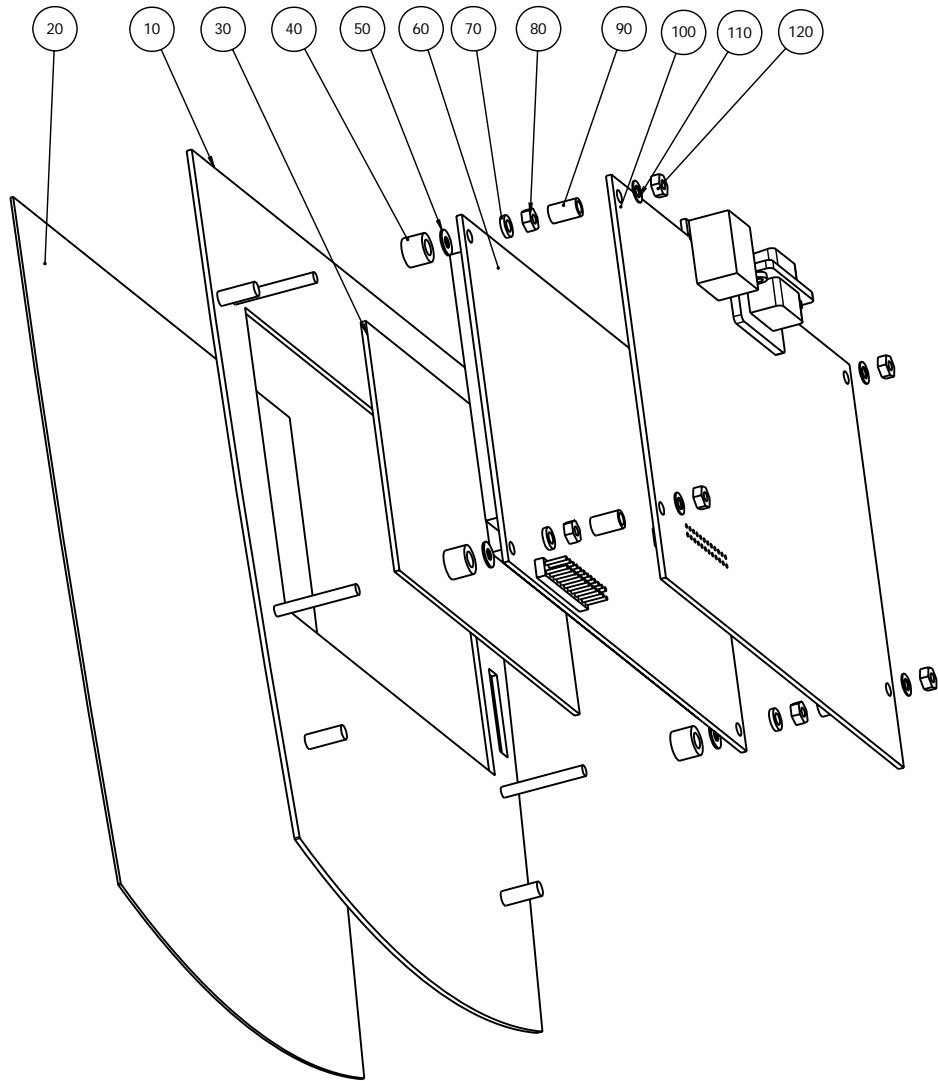
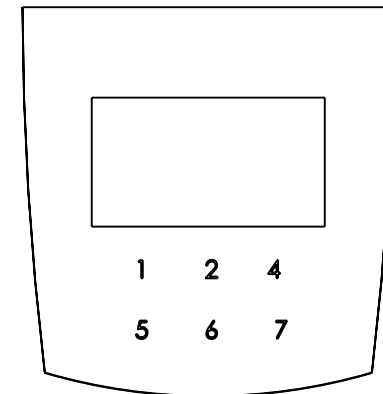
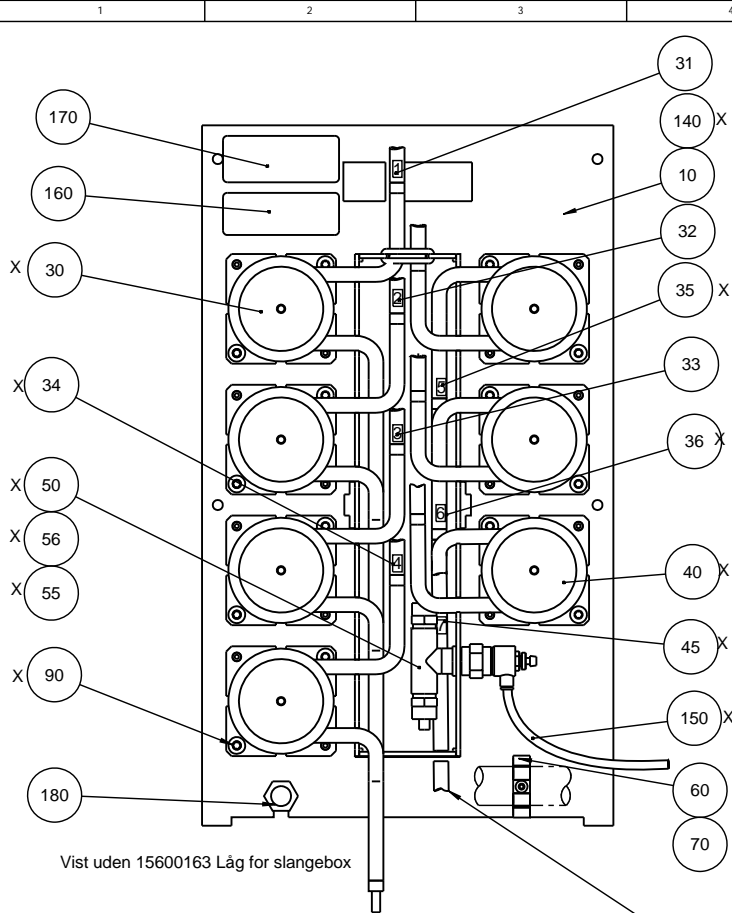
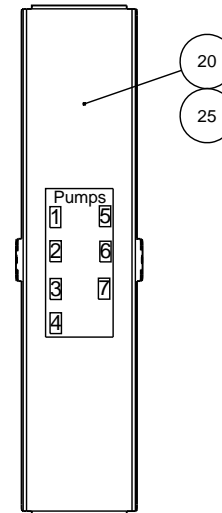
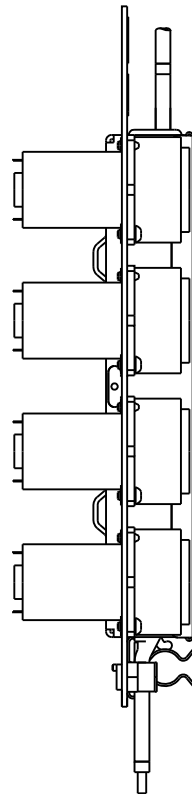
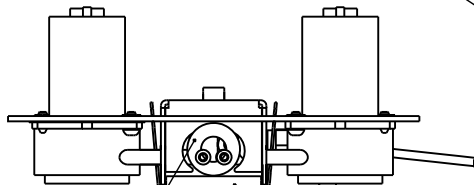
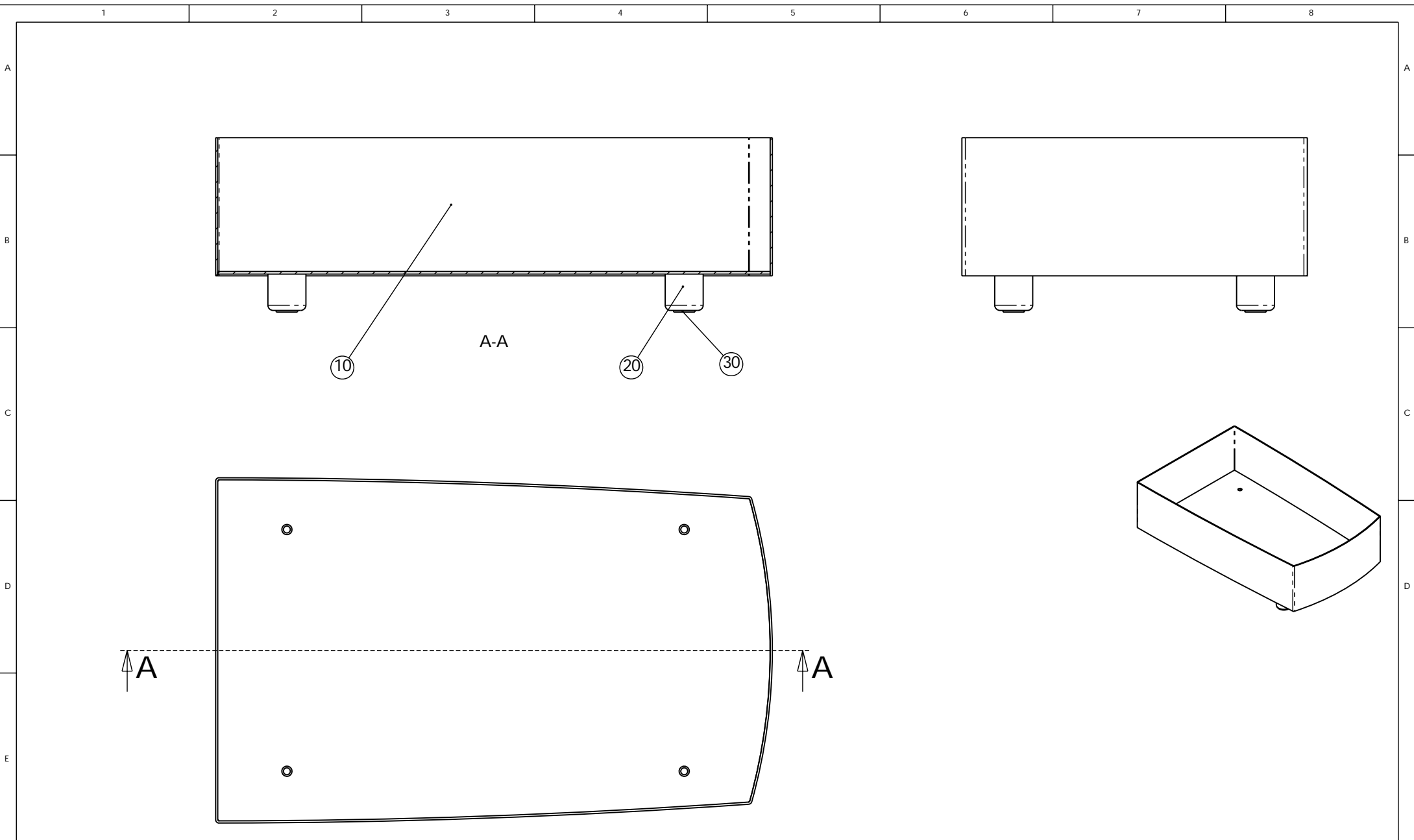
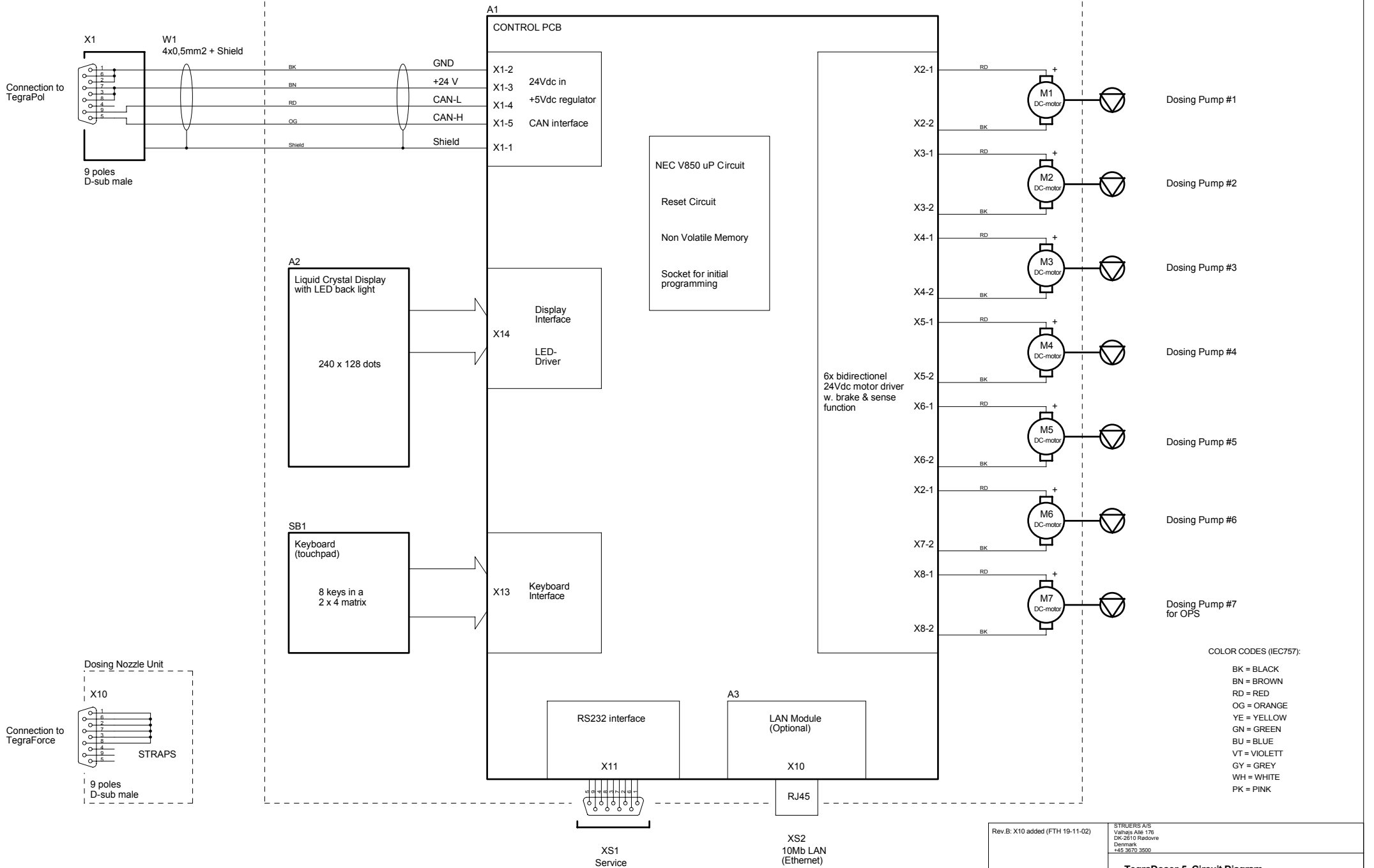
is in conformity with the provisions of the following directives:

Safety of Machinery 2006/42/EC according to the following standard(s):
EN ISO 12100-1:2005, EN ISO 12100-2:2005, EN ISO 14121-1:2007, EN 60204-1:2006/A1:2009,
EN 61010:2010.

EMC-Directive 2004/108/EC according to the following standard(s):
EN 61000-6-4:2007, EN 61000-6-1:2007, EN 61000-6-2:2005, EN 61000-6-3:2007, EN 61000-6-4:2002,
EN 61326-1:2006.

Low Voltage Directive 2006/95/EC according to the following standard(s):
EN 60204-1:2006/A1:2009, EN 61010:2010.

Supplementary Information The equipment also complies with the following standards:
UL3101-1:1993, NFPA70:2011, CAN-CSA 22.2 No. 1010-010-30,
FCC 47 CFR Part 15 Class A, AS/NZS 2064.1/2.

The above has been declared according to the global method, module A

Date: 31.01.2011

Christian Skjold Heyde,
Vice President, R&D and Production, Struers A/S

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Date: 31.01.2011

Christian Skjold Heyde,
Vice- President, R& D et Production, Struers A/S

The TegraSystem is designed as a flexible modular system, which you can combine to fulfil your specific needs. Use a single element or several. Expand and adapt the modules to meet your developing needs.

With the TegraSystem you will gain the benefits of efficiency, reproducibility and network integration.

And you will enjoy the reliability of proven technology from the market leader.

High efficiency and reproducibility is ensured due to precise control of all preparation parameters

Consumable costs are reduced and specimen throughput is increased due to improved efficiency

Easy to operate with large graphical display and pre-installed preparation methods

LED light illuminates the preparation area

The emergency switch stops all moving parts in case of an emergency

Networking allows you to share your methods and improve efficiency and standardization

Preparation methods are automatically recalculated to suit different system sizes

Only one power, water and air connection is necessary to supply the whole system

Durable and proven technology ensures years of trouble-free operation

The modular and flexible system

The Tegra family

With the TegraSystem you can choose between 2 grinding and polishing machines, 1 specimen mover and 2 dosing systems. The machines can be combined in several ways to suit your specific needs.

Reproducibility is the key

In all materialographic specimen preparation, being able to reproduce preparation results that are exactly the same every time is crucial to success. Many factors have an influence on your ability to achieve the necessary reproducibility, and many of these can be controlled by employing automated preparation equipment like the TegraSystem.

Ensuring reproducibility and efficiency

Your workload and need for flexibility should determine the size of your TegraSystem. However, no matter how many or few specimens you need to prepare, the investment in a TegraForce specimen mover will be beneficial. The operator avoids the often tedious manual preparation tasks, and since the TegraForce specimen mover controls the force applied precisely, specimens are prepared exactly the same way every time, independent of operator skills. For labs with many specimens, specimen movers are the only way to ensure cost-effectiveness since several specimens can be prepared simultaneously.

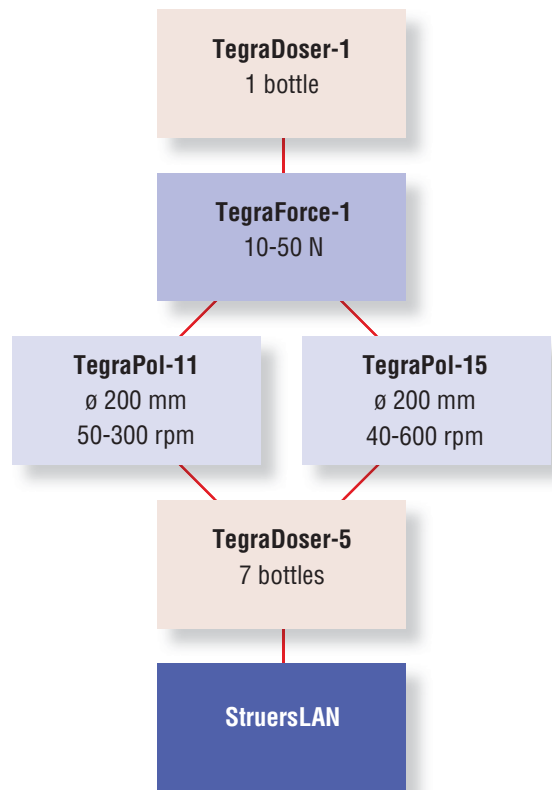
Adding a TegraDoser-5 dosing system further improves reproducibility, efficiency and cost-effectiveness. With the ability to store and recall preparation methods, this dosing system yields exactly the same preparation results every time the preparation procedure is performed.

The high level of automation also allows for unattended preparation, thus improving the lab efficiency by giving the operator free time for other tasks. Cost-effectiveness is further increased by automating and controlling consumable consumption.

The innovative approach

The modular concept lets your TegraSystem grow and develop with your needs and challenges. For the lab with few specimens or repetitive preparation of the same material every day, the TegraDoser-1 is an innovative approach to automated preparation.

TegraDoser-1 is a single dosing unit, which is attached directly onto the TegraForce specimen mover. The dosing unit not only controls the dosing of diamond suspension or lubricant precisely, but also controls the preparation parameters for that specific preparation step.

With 2-3 TegraDoser-1 units it is possible to run an entire preparation method simply by exchanging the dosing unit. Reproducibility is ensured since the preparation parameters are exactly the same every time. The necessary investment is low in relation to the benefits achieved.

Combined with the all-in-one diamond suspensions from Struers, this solution allows for unattended preparation and ensures high quality and cost-effective preparation results.

As you will see on the following pages, the TegraSystem holds unique advantages for fast, efficient and cost-effective grinding and polishing in a safe and healthy working environment.

The TegraPol grinding and polishing machines

The TegraPol grinding and polishing machines are a series of sturdy, reliable and powerful machines, based on proven technology and designed for preparation of all types of materialographic specimens.

- Can be used for manual preparation or be equipped with specimen mover and dosing system.
- Available in 2 speed versions for standard specimen preparation and special applications
- Easy and straightforward operation
- TegraPol automatically detects connected units, no configuration is necessary
- TegraPol power supplies all connected units
- Sturdy and powerful design ensuring years of trouble-free operation

The 200 mm TegraPol-11 and -15

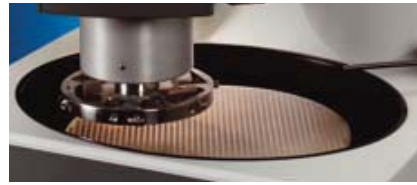
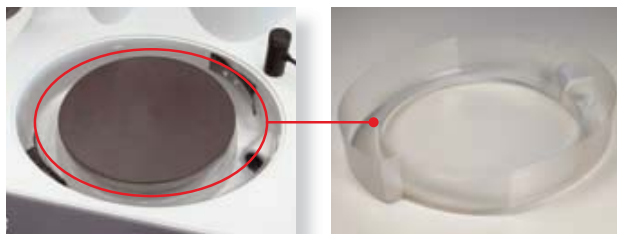
If you are preparing a limited number of specimens and still have high requirements to preparation quality, the 200 mm single disc machines are the perfect solution. TegraPol-11 is a unit covering all basic requirements with a speed of 50-300 rpm, whereas TegraPol-15 has a higher speed range from 40-600 rpm.

Added equipment is automatically detected

As soon as a TegraForce or a TegraDoser is connected, the TegraPol automatically detects which unit has been connected and no configuration is required. The operation of the whole system is carried out on the TegraPol. A large graphical display ensures easy and straightforward operation.

In the TegraSystem, the TegraPol is the only machine connected to the power mains. All other units are supplied directly from the TegraPol. The advantage here is of course the limited need for power outlets at your facilities. The power supply for all connected equipment is only 24 V.

For TegraPol-11 / -15 a disposable bowl liner is available to facilitate cleaning of the machine considerably.

The oval opening on the TegraPol machines grants easy access when exchanging consumables. At the same time it allows to place the specimen holder or specimen mover disc in the best possible position over the preparation disc to achieve optimum preparation results

Emergency stop

All TegraPol machines are equipped with an emergency stop, for immediate stopping of all moving parts in case of an emergency.

Environmentally-friendly recirculation system

A recirculation cooling unit can be connected to the TegraPol grinding/polishing machines. Controlled from the TegraPol, it is started and stopped automatically. Polluted debris from the preparation can be collected safely and the cooling water is utilised repeatedly. This reduces pollution and waste of water and eliminates the risk of corrosion of specimens by ordinary tap water.

Recirculation System

The TegraForce specimen mover

Mounting a specimen mover on the TegraPol is the first step towards automatic preparation. The precise control of all preparation parameters including the force is necessary to ensure reproducible preparation results. This is only obtained by electronic force control and, on TegraForce; the force on the specimens can be exerted very precisely, and the advantages are many:

- High reproducibility through precise control of preparation parameters
- Preparation of single specimens with individual force applied – flexibility is maintained and there is no need for dummy specimens

Preparation of single specimens

When single specimens are inserted in a specimen mover plate, the force is applied individually to each specimen. With this outstanding feature there is no need for dummy specimens. If your workflow demands short response-times you will benefit greatly from the individually applied force and maintain maximum flexibility. The preparation of individual specimens starts with grinding using fine grain sizes. There is no need for coarse pre-grinding since each specimen is levelled and ground independently. Very low force (10 N) can be applied on the individual specimen allowing for preparation of very fragile materials.

TegraForce-1

TegraForce-1 is designed to be mounted on a TegraPol-11 or TegraPol-15 and is suitable for preparation of both single specimens and specimens clamped in a specimen holder. 1, 2 or 3 specimens of up to 40 mm/1½" diameter can be fine ground and polished individually. For high removal rate or when plane grinding is required, the 120 mm diameter specimen holder where 3 specimens are prepared simultaneously is used.

To avoid damage when the preparation is started, the specimen holder or the pressure feet are lowered at reduced force. After a few seconds, the force applied to the specimens is slowly raised to the pre-set value.

On TegraForce the force can be reduced automatically towards the end of the preparation, resulting in lower material removal and smaller scratches. This saves time on the following preparation steps, since the amount of deformation to be removed is minimised.

TegraForce is equipped with LED illumination. While providing extra light to the preparation area, LED's grant very long life-time and minimum power consumption, making them most environmentally-friendly.

TegraForce-1

Accessories

Struers offers a wide range of accessories to the TegraForce-1 specimen mover. For TegraForce-1 both specimen mover plates and specimen holders are available, enabling the preparation of many different sizes and shapes of specimens. For the manufacture of your own special specimen holder designs, Struers can deliver specimen mover plates without holes.

LED illumination

The TegraDoser-5 dosing system

Adding an automatic dosing system to your TegraSystem gives you complete control of the entire preparation process and thus the quality of your preparation results and the costs involved.

- Built-in database with 10 preprogrammed Metalog Guide methods
- Built-in database with up to 200 methods leaves you in complete control
- Easy programming and large graphic display
- Excellent reproducibility
- Cost saving specimen preparation
- High-quality peristaltic pumps guarantee exactly the same dosing every time
- No vaporization or spray mist
- Automatic recalculation of method parameters when exchanged between different machines
- Automatic dosing of Oxide Polishing Suspensions
- Configuration according to your needs

Control the quality of your results

Manual dosing always implies certain variations in the amounts of abrasives and lubricants being dosed, and in the intervals between each dosing. TegraDoser-5 ensures that exactly the same amount is dosed each time – no matter how many times a preparation sequence is carried out.

But TegraDoser-5 is much more than just a dosing system. The built-in database can contain up to 200 preparation methods that can be stored and recalled at any time. This gives you full control of all preparation parameters and methods, regardless of which technician, machine or laboratory performs the job. Exactly repeatable settings give reliably reproducible results and frees resources from routine checks and repeats.

On a large graphical display preparation methods can easily be programmed and edited, whereas the method is selected and started from the TegraPol.

No vaporization or spray mist

The peristaltic pumps on the TegraDoser transport both suspensions and lubricants to the polishing cloth. In other dosing systems the suspension is vaporized and sprayed onto the cloth, resulting in an unhealthy spray mist. With the TegraDoser the liquid is dosed exactly where you want it on the cloth, improving preparation quality and improving the working environment.

Control your consumables costs

With the help of TegraDoser the dosing of suspensions and lubricants can be adjusted to the exact amount necessary for the preparation method in question. This will result in substantial savings on the consumables account. And TegraDoser also enables the operator to attend to other tasks during the preparation process – the dosing is automatically taken care of.

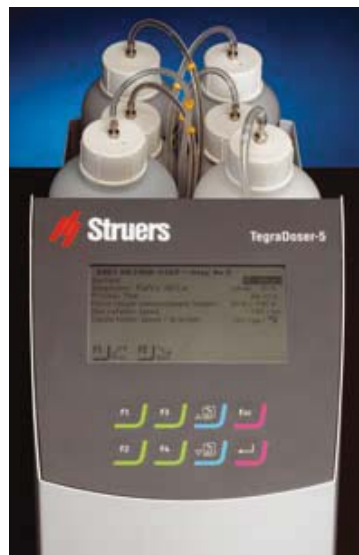
Automatic recalculation of preparation methods

A preparation method always needs to be adapted to the disc size of the grinding and polishing machine employed. As an example, the amount of suspension and lubricant that has to be applied on a 200 mm disc system is much less than the amount needed for a 300 mm disc system, simply because the surface area is smaller. The TegraSystem takes complete care of these required adjustments. Since the TegraSystem automatically detects the connected machines and units it can also automatically recalculate preparation methods to fit all system sizes. This ensures that preparation quality and reproducibility is maintained when preparation methods are being exchanged between machines, colleagues and labs.


Automatic dosing of Oxide Polishing Suspensions

The use of peristaltic pumps also makes it possible to include the dosing of oxide polishing suspensions in an automatic dosing unit. No other separate dosing units offer this feature. In addition, when working with oxide polishing suspensions, water is automatically applied onto the polishing cloth before and after the preparation. This ensures that the cloth is prepared prior to the polishing and that both specimens and cloth are cleaned afterwards.

TegraDoser-5 can be equipped with up to 7 peristaltic pumps, 6 for diamond suspensions or lubricants and 1 for Oxide Polishing Suspensions.

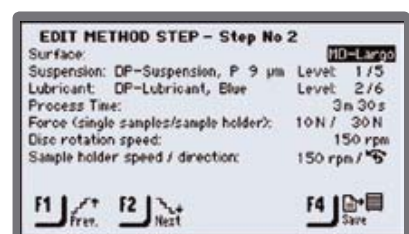Choose between 3, 4, 5 or 6 pumps for diamond suspensions/lubricant both with and without the pump for Oxide Polishing Suspensions

Data management and method details on TegraDoser-5

The TegraDoser-1 dosing system

If your specimen preparation only involves one method or a small number of specimens, the TegraDoser-1 concept can increase the automation level drastically at a low investment cost.

- High reproducibility and cost-efficiency are ensured due to precise dosing
- Individual ID settings let you dedicate one dosing unit for each preparation step
- Entire preparation methods can be run just by exchanging dosing units
- Especially well-suited for Struers all-in-one diamond suspension products allowing unattended specimen preparation

The single dosing unit TegraDoser-1 is attached directly onto the TegraForce and its built-in peristaltic pump enables precise control of the applied amount of suspension – your guarantee for constant dosing levels and high reproducibility.

The dosing levels are equivalent to the dosing levels of TegraDoser-5 and can be taken over directly. With the peristaltic pump the dosing levels are absolutely reproducible and all types of lubricants and diamond suspensions can be used.

Run an entire method with no programming at all

TegraDoser-1 is equipped with a small switch, which allows 10 different ID settings. The TegraPol recognizes the ID setting of each TegraDoser-1 and automatically recalls all the parameter settings related to the specific preparation step, which was run the last time the TegraDoser-1 in question was used. In this way, you can dedicate a TegraDoser-1 unit to each preparation step, and an entire preparation method can be run without any programming, simply by exchanging TegraDoser-1 units. This not only facilitates operation it also provides higher security because of increased reproducibility.

The winning combination of DiaPro and TegraDoser-1 lets you control the entire preparation method just by changing the dosing unit

Innovative all-in-one diamond products

Struers has developed innovative all-in-one DiaPro and DiaDuo diamond suspensions which combine diamond suspension and cooling lubricant into one product. This not only facilitates easy-of-use but also ensures that the consumables costs are kept under control. These products are extremely well-suited for the TegraDoser-1 dosing system as they will allow for unattended specimen preparation.

For more information about Struers diamond products, please ask for separate product brochures.

The StruersLAN networking module

The TegraDoser-5 can be equipped with a StruersLAN module. This allows you to connect your TegraSystem to your company network or to a single PC.

- Easier programming of preparation methods from any connected PC in the company
- Preparation methods can be exchanged between machines, colleagues and labs
- Standardisation and efficiency is improved instantly
- Preparation methods and consumables lists can be downloaded from the Internet and used immediately

Networking simplifies programming

The StruersLAN module is a network card, which works according to the Ethernet standard, used worldwide in data communication. By providing the TegraDoser-5 with the StruersLAN module, the TegraSystem can be connected to a single PC or directly to the company's LAN network (Local Area Network), allowing the TegraSystem to be accessed from all the connected PCs.

In both cases programming of preparation methods is made extremely easy, because a complete keyboard is available. Also copy and paste functions, rename, delete, insert, - in short - all the functions usually available on a PC are available here too. The program is based on the same functions as the Windows™ Explorer.

Networking lets you share your knowledge

Do you want to share your preparation method with a colleague, another lab in your company or even another company? With the TegraSystem this is possible.

All TegraSystem modules connected to the same network can be accessed by PCs connected, and preparation methods can easily be copied directly from one machine to another.

Differences in machine size and set-up are not a problem since the TegraSystem automatically recalculates all preparation parameters when a method is transferred from one system to another.

Even when the TegraSystems are not connected to the same network, preparation methods can be exchanged. Simply export the preparation method, send it attached to an e-mail and import it into another TegraSystem without any problems.

Networking improves standardisation and efficiency

The networking possibilities greatly improve visibility and efficiency. Preparation methods can easily be made accessible on all the TegraSystems in the lab. Redundant work performed in different labs within the company can be eliminated. A central lab can standardise company methods and easily transfer them to TegraSystems in all satellite labs, ensuring consistent settings throughout the organisation. As soon as a method is approved, it can be transferred and made accessible instantly. Even specimen preparation performed at sub-suppliers can be standardised. The possibilities are numerous as required to match the different needs in modern companies.

Download of methods from the Internet

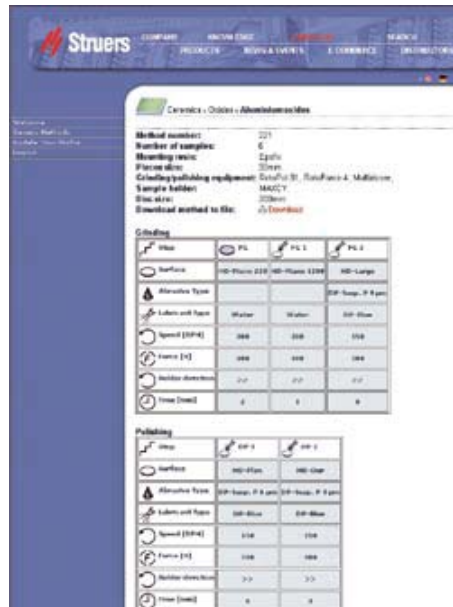
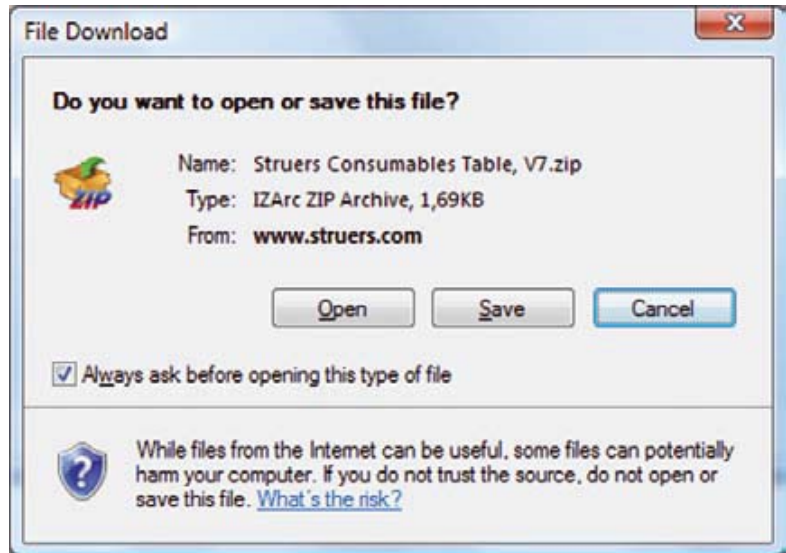
With a TegraSystem and StruersLAN networking possibilities, you benefit from the Internet tools provided by Struers. Hundreds of different preparation methods can be downloaded directly from Struers e-Metalog on the Internet in a format that is compatible with the TegraSystem. The methods can be employed instantly or adjusted to fit your specific application.

Also available from our homepage is the consumables table used by the TegraSystem. Whenever we release new grinding discs, polishing cloths, diamond suspension etc. the list will be updated and is immediately available for download. The advantage of this list lies in the built-in interaction between the consumables. If you e.g. select SiC as abrasive, the lubricant is automatically set to water and it will not be possible to select a suspension or a lubricant.

Full network integration

State-of-the-art communications technology gives you the potential to synchronise machine settings and preparation methods throughout your laboratory network. You can implement changes the moment they are approved. Using a central database simplifies the standardisation processes and keeps you in full control.

The flexibility of the TegraSystem allows you to tailor the system to meet your exact demands. Use a single element or several. Grow and adapt the modules to meet your developing needs. Gain the benefits of reproducibility and efficiency – in other words, let the TegraSystem revolutionize your world of sample preparation.

Hundreds of different preparation methods can be downloaded directly from Struers e-Metalog on the Internet.

State-of-the-art communications technology gives you the potential to synchronise machine settings and preparation methods throughout your laboratory network.

Grinding/Polishing machine	TegraPol-15	TegraPol-11
Disc diameter	200 mm / 8"	
Disc speed	40 - 600 rpm	50 - 300 rpm
Rotational direction	CCW	
Power supply	1/3 200-240 V / 50-60 Hz	
Motor Power continuous / max	370 / 555 W	180 / 270 W
Automatic torque control	Yes	
Emergency stop	Yes	
Touch Pad control	Yes	
Graphical display for easy selection of preparation methods	Yes	
Soft start with pre-dosing	Yes	
Controlled acceleration and deceleration	Yes	
High speed drying of preparation disc, Spin Function @ 600 rpm	Yes	
Max. torque at disc	>20 Nm	>10 Nm
Compatible with Struers MD-system	Yes	
Automatic detection of connected equipment	Yes	
Connection of a polishing head	TegraForce-1	
Maintenance management	Diagnostic and time used	
Service features	Machine test and self diagnostic	
Width, Depth, Height and Weight	410 mm / 670 mm / 315 mm / 25 kg	410 mm / 670 mm / 315 mm / 23 kg

Polishing head	TegraForce-1
Single specimens	Yes, 1 - 3 specimens
Specimen holder	Yes, 3 specimens
Rotational speed	50 - 150 rpm
Rotational direction	CCW + CW
Motor Power continuous / max (power supplied from TegraPol)	35 / 53 W
Force per specimen, single specimens / specimen holder	10 - 50 N / 30 - 150
LED lighting	Yes
Lifting capacity for specimen holders	-
Connection of dosing system	TegraDoser-5 or TegraDoser-1
Identical Start and Stop position of polishing head	Yes
Soft start function	Yes
Force reduction	Yes
Automatic lift-off after preparation	No
Locking of polishing head	Manual
Push button for quick release of specimen holders	No
Guide for positioning of polishing head	Yes
Fixed stop position of polishing head	Yes
Service features	Machine test and self diagnostic
Width, Depth, Height and Weight	140 mm / 311 mm / 445 mm / 7 kg

Dosing unit	TegraDoser-5	TegraDoser-1
Precise dosing with peristaltic pumps	Yes	
Number of pumps	Up to 7 pumps	Single pump
Handles lubricant, suspension and all-in-one suspensions	Yes	
Handles OP-suspension	Yes	-
Integrated automatic cleaning function	Yes	-
Graphical display for programming of preparation methods	Yes	-
Integrated database	10 Metalog Guide methods and 200 user methods	-
ID switch for parameter step control	-	Yes
Download of preparation methods from Internet (Struers e-Metalog)	Yes	-
Possibility to program methods via PC	Yes	-
Access to consumables list	Yes	-
Width, Depth, Height and Weight	200 mm / 210 mm / 380 mm / 8.5 kg (excl. bottle tray)	120 mm / 50 mm / 240 mm / 550 g

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CANADA

Struers Ltd.
7275 West Credit Avenue
Mississauga, Ontario L5N 5M9
Phone +1 905-814-8855
Fax +1 905-814-1440
info@struers.com

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Struers' equipment is in conformity with the provisions of the applicable International Directives and their appurtenant Standards. (Please contact your local supplier for details)

Struers' products are subject to constant product development. Therefore, we reserve the right to introduce changes in our products without notice.

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All rights reserved. © Struers 20F€
Struers Inc.
Struers GmbH